



Universiteit
Leiden
The Netherlands

Development of a kidney-on-a-chip model for compound screening and transport studies

Vormann, M.K.

Citation

Vormann, M. K. (2021, September 9). *Development of a kidney-on-a-chip model for compound screening and transport studies*. Retrieved from <https://hdl.handle.net/1887/3209238>

Version: Publisher's Version

License: [Licence agreement concerning inclusion of doctoral thesis in the Institutional Repository of the University of Leiden](#)

Downloaded from: <https://hdl.handle.net/1887/3209238>

Note: To cite this publication please use the final published version (if applicable).

Cover Page



Universiteit Leiden



The handle <https://hdl.handle.net/1887/3209238> holds various files of this Leiden University dissertation.

Author: Vormann, M.K.

Title: Development of a kidney-on-a-chip model for compound screening and transport studies

Issue Date: 2021-09-09



Chapter 3

Implementation of a Human Renal Proximal Tubule on a Chip for Nephrotoxicity and Drug Interaction Studies

MK Vormann*, J Vriend*, HL Lanz, L Gijzen, A van den Heuvel, S Hutter, J Joore, SJ Trietsch, C Stuut, TTG Nieskens, JGP Peters, D Ramp, M Caj, FGM Russel, B Jacobsen, A Roth, S Lu, JW Polli, AA Naidoo, P Vulto, R Masereeuw, MJ Wilmer, L Suter-Dick

*Authors contributed equally

Based on:

Journal of Pharmaceutical Sciences 110, 1601-1614. (2018)

Abstract

Proximal tubule epithelial cells (PTEC) are susceptible to drug-induced kidney injury (DIKI). Cell-based, two-dimensional (2D) *in vitro* PTEC models are often poor predictors of DIKI, probably due to the lack of physiological architecture and flow. Here, we assessed a high throughput, 3D microfluidic platform (Nephroscreen) for the detection of DIKI in pharmaceutical development. This system was established with four model nephrotoxic drugs (cisplatin, tenofovir, tobramycin and cyclosporin A) and tested with eight pharmaceutical compounds. Measured parameters included cell viability, release of lactate dehydrogenase (LDH) and N-acetyl- β -D-glucosaminidase (NAG), barrier integrity, release of specific miRNAs, and gene expression of toxicity markers. Drug-transporter interactions for P-gp and MRP2/4 were also determined. The most predictive readouts for DIKI were a combination of cell viability, LDH and miRNA release. In conclusion, Nephroscreen detected DIKI in a robust manner, is compatible with automated pipetting, proved to be amenable to long-term experiments, and was easily transferred between laboratories. This proof-of-concept-study demonstrated the usability and reproducibility of Nephroscreen for the detection of DIKI and drug-transporter interactions. Nephroscreen it represents a valuable tool towards replacing animal testing and supporting the 3Rs (Reduce, Refine and Replace animal experimentation).

Introduction

Renal proximal tubules are susceptible to drug-induced kidney injury (DIKI) [1], which can be a dose-limiting factor in pharmacotherapy and lead to kidney failure in patients. Up to 14–26% of acute kidney injury cases are caused by DIKI [2], [3], [4]. This clinical observation is in stark contrast to the low number of drug candidates, 2% of total, that fail in preclinical development due to nephrotoxicity [5]. This disparity demonstrates the need for better clinically predictive models for nephrotoxicity for use in discovery and early development. Screening-friendly, human-relevant *in vitro* test systems will also help replace animal studies in drug discovery and contribute to the 3 Rs (Replacement, Reduction and Refinement of animal studies). One third of drugs and drug candidates tested are (partially) excreted via the urine following dosing to humans [6], [7]. Active secretion of xenobiotics in the kidney takes place mainly via the proximal tubule epithelial cell (PTEC) of the nephron [8]. A polarized monolayer of PTECs, joined by characteristic epithelial tight junctions, separates the vasculature on the basolateral side from the tubular fluid on the apical side and regulates solutes and water by active transport mechanisms.

In recent years, three-dimensional (3D) microfluidic *in vitro* models of PTECs, also referred to as proximal tubule-on-a-chip, gained significant interest as predictive platforms for nephrotoxicity in drug development [9]. Currently employed *in vitro* two-dimensional (2D) PTEC models lack important *in vivo* characteristics, such as cell-extracellular matrix (ECM) interaction and fluid shear stress (FSS), which limits their relevance and predictivity with regards to nephrotoxicity [1],

[9], [10]. Renal proximal tubule-on-a-chip showed improved characteristics, such as increased tight-junction formation (ZO-1 expression), and increased number of cilia and microvilli at the apical membrane [11], [12]. Features like albumin uptake and increased P-glycoprotein (P-gp) activity in renal proximal tubule-on-a-chip, and nephrotoxicity induced by cisplatin showed close resemblance with *in vivo* observations [11], [12].

Implementation of renal proximal tubules-on-a-chip in large-scale nephrotoxicity screening is limited due to the complexity and low throughput of most models, often consisting of one chip connected to pumps to generate flow [11], [12]. Furthermore, most renal proximal tubules-on-a-chip described lack basolateral and apical compartments or make use of a two-compartmental model separated by cells cultured on a semi-permeable membrane ignoring cell-ECM interaction [12]. Choice of renal cell source is another important factor in the proximal tubule-on-a-chip. Although freshly isolated primary PTECs show more physiological characteristics in a nephrotoxicity screening proximal tubule-on-a-chip model [13], availability of primary PTECs limits the throughput of this model. Using immortalized renal PTECs would not only overcome this problem but would also enhance reproducibility across different laboratories [13].

This study aimed to assess a high throughput, 3D-microfluidic platform (Nephroscreen) for the detection of drug-induced nephrotoxicity: This platform was specifically designed to fulfill requirements of pharmaceutical companies and promote alternative methods to animal testing in support of the 3Rs. Nephroscreen was first challenged with four selected model human nephrotoxic drugs (cisplatin, tenofovir, tobramycin and cyclosporin A) that affect the proximal tubule. This was complemented with eight additional compounds provided by three pharmaceutical companies for evaluation in a blinded manner. The data were generated at several laboratories to ensure that the platform is robust and transferable.

In Nephroscreen we combined a microfluidics platform with suitable cell lines and appropriate assays. An automatable, microfluidic plate consisting of multiple chips, the OrganoPlate, was combined with renal PTEC cell lines, exposed to FSS induced through passive levelling by gravity [14], [15], [16]. The chosen cells, conditionally immortalized PTEC overexpressing OAT1 (ciPTEC-OAT1) or pseudo-immortalized renal PTEC (RPTEC), are able to establish a polarized epithelium expressing functional transporters [15], [16]. The ciPTEC-OAT1 has been engineered to overexpress OAT1, in addition to other transporters, such as P-gp and OCT2, expressed in the parental ciPTEC line [17], [18]. Therefore, they are an ideal cell line to study the toxicity of substances requiring transport into the cell via these transporters. RPTEC, on the other hand, were chosen due to their performance in establishing leak-tight epithelial barriers. Thus, both implemented cell lines show different, complementary strengths that were exploited by measuring suitable endpoints: ciPTECs were utilized for cytotoxicity, biomarker, and transporter activity measurements while RPTECs performed very well in functional assays such as barrier function.

Cellular damage was studied using various read-outs, such as enzymatic production of formazan (WST-8 assay), indicator of cell viability and release of lactate dehydrogenase (LDH) and N-acetyl- β -D-glucosaminidase (NAG), measures for membrane integrity. Molecular markers of cellular stress included extracellular levels of specific miRNAs: mir34a, mir21, mir192, and mir29a [14], [15], as well as gene expression of heme oxygenase 1: HMOX1 and neutrophil gelatinase-associated lipocalin: NGAL [19]. Functional parameters for the epithelial monolayer consisted of barrier integrity testing as well as interactions with P-gp and multidrug resistance-associated proteins 2 and 4 (MRP2/4) [15], [16].

This study was part of the NC3Rs (<https://www.nc3rs.org.uk/>) crackIT challenge Nephrotube, launched to generate predictive, animal free-systems for the detection of nephrotoxicity. The results showed the potential of the Nephroscreen for predictive, animal free-detection of nephrotoxicity and drug-transporter interactions.

Materials and Methods

Experimental Workflow: Multi-laboratory Collaboration

The combination of assays described previously [14], [15], [16], used towards a functional Nephroscreen, were performed in laboratories at the School of Life Sciences, University of Applied Sciences Northwestern Switzerland in Muttenz, Switzerland (lab B), the department of Pharmacology and Toxicology, Radboud University Medical Center, Nijmegen, The Netherlands (lab N) and Mimetas, Leiden, The Netherlands (lab L). Transferability, robustness and reproducibility of Nephroscreen was established by performing experimental procedures in the three laboratories using two renal cell lines (ciPTEC and RPTEC). The experimental workflow is depicted in Fig. 1. For cell viability assessment, tubules of both cell lines were cultured and then exposed to nephrotoxics for 24 or 48 h (lab L). Measurement of cell viability included enzymatic production of formazan using WST-8 (Lab L), LDH release (lab N) and gene expression of toxicity markers (lab L and N). In addition, RPTEC were tested for their barrier function (Lab L). In lab N, the effect of the nephrotoxics on the transporter functionality of P-gp and MRP2/4 was assessed in ciPTEC-OAT1. The release of selected miRNAs (mir-21, -29a, -34a and -192) and secretion of the enzyme N-acetyl- β -D-glucosaminidase (NAG) in ciPTEC-OAT1 were measured in lab B.

Cell Culture

CiPTEC-OAT1 cells (Cell4Pharma, Oss, NL) were cultured in T75 flasks in a 1:1 mixture of Dulbecco's modified Eagle's medium and nutrient mixture F-12 without phenol red (DMEM-F12, Gibco, ThermoFisher Scientific, Carlsbad, USA, 11039-021) supplemented with insulin-transferrin-sodium selenite media supplement (Sigma-Aldrich, I1884, insulin 5 μ g/mL; transferrin 5 μ g/mL; sodium selenite 5 ng/mL), 36 ng/mL hydrocortisone (Sigma-Aldrich, H0135), 10 ng/mL epidermal growth factor (Sigma-Aldrich, E9644), 40 pg/mL 3-iodothyronine (Sigma-Aldrich, T5516),

10% v/v fetal bovine serum (FBS, Gibco, 16140-071 or Greiner Bio-One, Alphen aan den Rijn, The Netherlands, 758093), and 1% v/v penicillin/streptomycin (P/S, Sigma-Aldrich, P4333), referred to as ciPTEC complete medium.

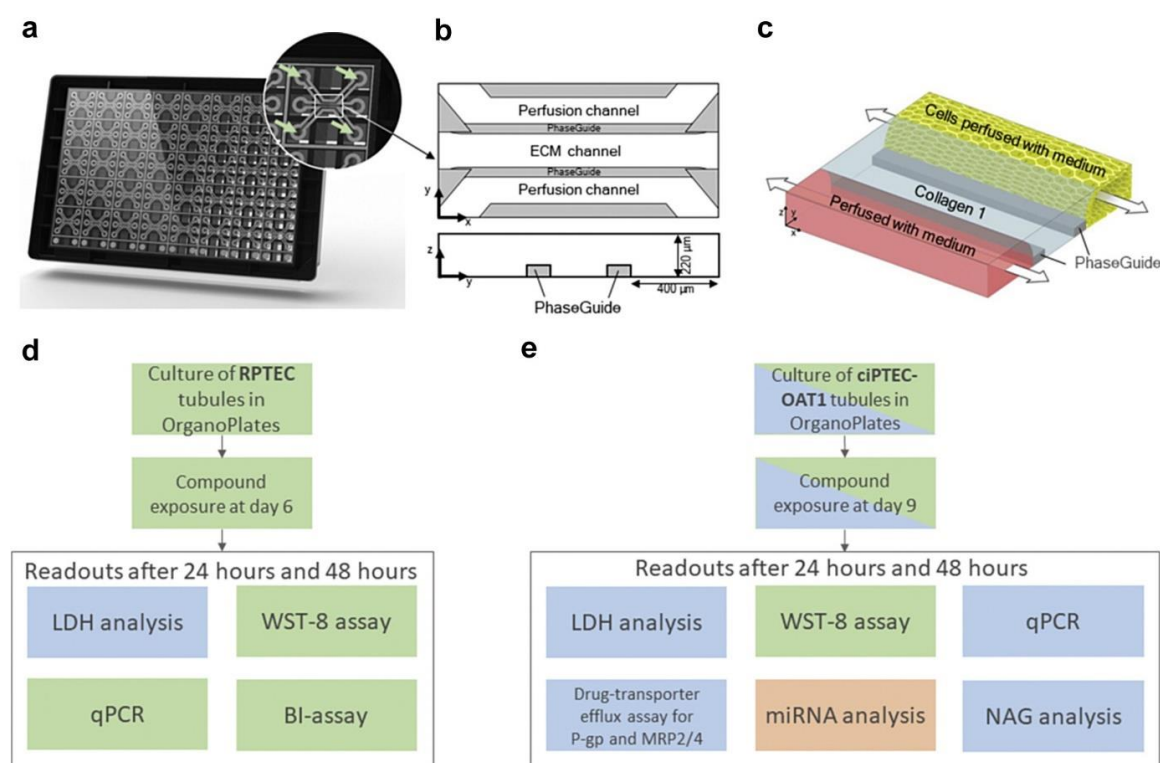


Figure 1: Overview of the microfluidics platform, the cell culture, exposure, and read-out methods performed on ciPTEC-OAT1 tubules and RPTEC tubules. (a) Image of the back side of the OrganoPlate 3-lane. The microfluidic network is positioned in-between a glass sandwich of two microscope grade glass plates which are attached to the bottom of a standard 384 titer well plate. Access to the microfluidic system is facilitated via the top wells. One OrganoPlate comprises in total 40 chips as 8 “wells” compose one chip. Green arrows indicate the inlets used for compound dosing. (b) Schematic of one chip presenting two perfusion channels and the extracellular matrix (ECM) channel in the middle. Channels are divided by 55 µm high ridges called PhaseGuide (grey bars) which act as pressure barriers. (c) Artist impression of one chip. The chip was loaded with collagen 1 (blue) to the ECM channel and proximal tubule cells (yellow) were seeded to the top channel. After cell attachment medium was added to both perfusion channels and perfusion was started (indicated by white arrows). (d) and (e) Flow charts indicating the culture of RPTEC and ciPTEC-OAT1, respectively. Each cell line was cultured following optimized conditions before exposure to compounds for 24 and 48 h. (d) Readouts for RPTEC included lactate dehydrogenase (LDH)-release, cell viability assay using WST-8 assay, determination of barrier integrity (BI) and collection of total RNA for quantitative PCR (qPCR). (e) Readouts for ciPTEC-OAT1 included lactate dehydrogenase (LDH)-release, N-acetyl-β-d-glucosaminidase (NAG) release, cell viability assay using WST-8 assay, determination of drug-transporter efflux assay, collection of release RNA for miRNA determination, and collection of total RNA for quantitative PCR (qPCR). Assays were performed and optimized in three laboratories and the obtained results are colour coded: Lab B, orange, Lab L (green), Lab N (blue).

Cells were incubated in a humidified incubator (33 °C, 5% v/v CO₂) and every 2–3 days ciPTEC complete medium was changed. At 90–100% confluency cells were washed with HBSS (Sigma-

Aldrich, H6648 or Gibco, 14025-100), detached with accutase (Sigma-Aldrich, A6964), pelleted (200–300×g, 5 min), and used for passaging (10,000–20,000 cells/cm²) or for seeding in the OrganoPlate. Cells were used for experiments between passage 52 and 65.

RPTEC cells (Kidney PTEC Control Cells, SA7K Clone, Sigma-Aldrich, Schnelldorf, Germany, MTOX1030) were seeded in T75 cell culture flasks coated with PureCol (Advanced BioMetrix, 5005-B, diluted 1:30 in cold Hank's balanced salt solution (HBSS, Sigma-Aldrich, H6648), 20-min incubation at 37 °C) in MEME alpha Modification (Sigma-Aldrich, M4526) supplemented with RPTEC Complete Supplement (Sigma-Aldrich, MTOXRCSUP), L-glutamine (1.87 mM, Sigma-Aldrich, G7513), Gentamicin (28 µg/mL, Sigma-Aldrich, G1397) and Amphotericin B (14 ng/mL, Sigma-Aldrich, A2942), referred to as RPTEC complete medium. Cells were incubated in a humidified incubator (37 °C, 5% v/v CO₂) and every 2–3 days RPTEC complete medium was changed. At 90–100% confluency cells were washed with HBSS, detached with accutase (Sigma-Aldrich, A6964), pelleted (140×g, 5 min), and used for seeding in the OrganoPlate. Cells for experiments were used up to passage 3.

OrganoPlate Culture

For all experiments, plates were seeded as described in our previous work [14, [15], [16]. Briefly, a three-lane OrganoPlate (Mimetas BV, 4003 400B, Fig. 1) with a channel width of 400 µm and a height of 220 µm was loaded with 1.6–2 µL of ECM gel composed of 4 mg/mL collagen I (AMSBio Cultrex 3D Collagen I Rat Tail, 3447-020-01), 100 mM HEPES (Life Technologies, ThermoFisher Scientific, 15630), and 3.7 mg/mL sodium bicarbonate (Sigma-Aldrich, 320 S5761) to the middle inlet of all 40 chips. After polymerization of the ECM, HBSS was added on top of the collagen I and the plate was incubated in a humidified incubator (37 °C, 5% v/v CO₂) overnight. RPTEC or ciPTEC-OAT1 were detached from culture flasks and resuspended at a concentration of 10×10^6 cells per mL in RPTEC complete medium or 20×10^6 cells per mL in ciPTEC complete medium, respectively. Of the cell suspension, 2 µL was injected into each top inlet, followed by an addition of 50 µL medium to the same well. For control chips, no cell suspension was added. Subsequently, the OrganoPlate was placed for 4 h at an angle of about 75° into a humidified incubator (37 °C, 5% v/v CO₂ for RPTEC or 33°, 5% v/v CO₂ for ciPTEC-OAT1). After attachment of the cells, which was confirmed by light microscopy, 50 µL medium was added to the top outlet, bottom inlet and bottom outlet. The OrganoPlate was placed flat in an incubator (37 °C, 5% v/v CO₂ for RPTEC or 33 °C, 5% v/v CO₂ for ciPTEC-OAT1) on an interval rocker platform (±7° angle, 8 min interval) enabling a bidirectional flow through the perfusion channels. At day 3, antibiotic free medium was used to refresh the medium in the chips, and the ciPTEC-OAT1 plates were transferred to a humidified incubator (37 °C, 5% v/v CO₂).

Compound Exposure

Cells were exposed for 24 or 48 h to four model nephrotoxics and eight blinded compounds following the general experimental design depicted in Fig. 1. The selection of the test concentrations was based on doses where toxicological effects were seen *in vivo* and on preliminary cytotoxicity data in 2D obtained with ciPTEC-OAT1. Briefly, cells were seeded in 96-well flat-bottom plates at 35000 cells/cm², maintained until confluency, subsequently exposed for 24 or 48 h to a concentration range of each of the compounds diluted in ciPTEC complete medium and viability was determined (WST-8 assay). For Nephroscreen testing, a low and a high toxicity concentrations were selected. In absence of detectable toxicity, the two highest concentrations were chosen (data not shown). Detailed solvent, dilutions, and culture conditions for each substance are described in Table S1, while data on previous studies used for the determination of potential test concentrations are listed in Table S2, together with relevant toxicity information disclosed after data acquisition and analysis of data had been completed. For compounds R1 and R2, additional experiments including more concentrations and additional time points were performed in lab B. For R1, cells were exposed for 48 h and for R2, cells were exposed for 11 days, with medium changes every 2–3 days.

Using a pipette, 50 µL of the model compounds (cisplatin, tenofovir, tobramycin and cyclosporin A) or the eight blinded compounds (G1, G2 and G3; R1, R2 and R3; and P1 and P2) diluted in medium were dosed via each of the 4 inlets guiding to the microfluidic channels (details in Fig. 1). Medium from both the apical and basolateral perfusion channels was refreshed before exposure to the nephrotoxics. For RPTEC, the exposure was started at day 6 and dilutions of nephrotoxics were prepared in TOX medium (MEME alpha Modification (Sigma-Aldrich, M4526) supplemented with RPTEC Tox Supplement (Sigma-Aldrich, MTOXRTSUP), L-glutamine (1.87 mM, Sigma-Aldrich, G7513). For ciPTEC-OAT1, cells were exposed on day 9 and dilutions of nephrotoxics were prepared in ciPTEC complete medium.

WST-8 Assay

Cell viability was determined using the Cell Counting Kit-8 (WST-8, Sigma-Aldrich, 96992) as described previously [14, [15]. Briefly, the WST-8 solution was diluted with ciPTEC complete medium or RPTEC TOX medium (ciPTEC-OAT1 or RPTEC respectively) and added to the channels of the OrganoPlate. After 20–60 min, absorbance was measured at 450 nm with a Multiskan™ FC Microplate Photometer (ThermoFisher Scientific) or with a Benchmark Plus microplate spectrophotometer (Bio-Rad, Veenendaal, The Netherlands). Viability was expressed as % of the control.

LDH Activity in Medium

LDH activity in the culture medium is a measure for membrane integrity. Medium was collected on ice after compound exposure and stored at –80 °C until further processing. After thawing, LDH activity was measured using an activity assay kit (Sigma-Aldrich, MAK066), following the

manufacturer's protocol. In brief, 5 μ L per sample was added to a flat bottom 96-well plate. In addition, a calibration curve using a NADH (1.25 mM) standard was prepared. Assay buffer was added to a final volume of 50 μ L per well and then a master reaction mix was added per well (1:1, v/v). After 3 min, absorbance was measured at 450 nm on a Benchmark Plus microplate spectrophotometer (Bio-Rad) every 5 min until absorbance measured in a sample was higher than the highest level of NADH in the calibration curve (12 nmol/well). Extracellular LDH activity was expressed as mU/mL.

NAG Measurement in Medium

The levels of β -N-acetylglucosaminidase (NAG) in supernatant are a measure for membrane integrity and are also used as a biomarker of kidney damage *in vivo*. They were determined using the NAG-Assay Kit (BioVision, K733-100), following provider's instructions. Briefly, 30 μ L of the supernatant (cell culture medium) were adjusted to a volume of 70 μ L with NAG assay buffer in a 96-well plate. To this, 55 μ L substrate were added and incubated for 30 min at 37 °C, followed by the addition of 25 μ L stop solution. After another 10 min at 37 °C, absorbance was measured at 400 nm in a Flexstation 3 (Molecular Devices). Absorbance of a standard dilution series of pNP (0–20 nmol of pNP/well) was measured in parallel and used for the calculations. Extracellular NAG activity was calculated (mU/mL) and expressed as % of the control.

Barrier Integrity Assay

The barrier integrity assay was performed as previously published on RPTEC [15]. Briefly, the medium of the perfusion channel (apical side of the cells) was replaced with RPTEC TOX medium containing 0.5 mg/mL tetramethylrhodamine (TRITC)-dextran (4.4 kDa, Sigma-Aldrich, FD20S) and 0.5 mg/mL fluorescein isothiocyanate (FITC)-dextran (150 kDa, Sigma-Aldrich, T1287). Subsequently, plates were imaged with the ImageXpress Micro XLS-C High Content Imaging System (Molecular Devices). Leakage of the dyes from the lumen (apical compartment) to the basal compartment into the ECM was measured over time, and the ratio between the basal and the apical compartment was analyzed with Image J version 1.51n [20]. From these measurements, the apparent permeability index (P_{app}) was calculated as previously described [15], using the following formula:

$$P_{app} = \frac{\Delta C_{receiver} \times V_{receiver}}{\Delta t \times A_{barrier} \times C_{donor}} \left(\frac{cm}{s} \right)$$

$\Delta C_{receiver}$ is the measured normalized intensity difference of the ECM to the donor channel (apical compartment) at t_{0min} and t_{10min} , $V_{receiver}$ is the volume of the measured region in the ECM channel (0.0001 cm³), Δt is the time difference (10 min), $A_{barrier}$ is the surface of the ECM interface with the medium channel (0.0057 cm²), and C_{donor} is the donor concentration of the dextran dyes (0.5 mg/mL).

Gene Expression of Toxicity and Nephrotoxicity Markers

Total intracellular RNA was isolated from cells harvested from perfusion channels in the OrganoPlate using the RNeasy Micro Kit (Qiagen, Venlo, The Netherlands) as previously described [16]. Complementary DNA (cDNA) was synthesized with Moloney Murine Leukemia Virus (M-MLV, Promega, M1705) reverse transcriptase, according to the manufacturer's protocol. Gene expression levels of heme oxygenase 1 (HMOX1) and neutrophil gelatinase-associated lipocalin (NGAL), encoded by the lipocalin-2 (LCN2), were measured using quantitative PCR (qPCR) with GAPDH as reference. TaqMan Universal PCR Master Mix (Life Technologies, 4304437) and gene specific primer-probe sets (HMOX1: Hs01110250_m1, LCN2: Hs01008571_m1 and GAPDH: Hs99999905_m1) were purchased from Applied Biosystems (ThermoFisher Scientific). The real-time PCR was carried out using a fluorometric thermal cycler (Qiagen, Rotor-Gene Q 3000 or Celtic Diagnostics, Corbett research PCR Rotorgene 6000). The PCR program consisted of 10 min of initial denaturation at 95 °C followed by 45 cycles of 15 s at 95 °C and 1 min at 60 °C. Fluorescence was detected at the end of each cycle at 510 nm (excitation 470 nm). Ct values were determined using the second derivative method. For each sample, dCt was calculated using Ct values of the gene of interest and the housekeeping gene (GAPDH). Treatment induced differential gene expression was calculated using the $-ddCt$ s $-(dCt_{\text{treatment}} - dCt_{\text{control}})$ and fold changes as 2^{-ddCt} .

Detection of miRNAs in Medium

Total RNA was extracted from 80 μ L supernatant using the miRNeasy® Serum/Plasma Kit (Qiagen, 217184) following the manufacturer's protocol. miRNAs were reverse-transcribed using the TaqMan® MicroRNA Reverse Transcription Kit (Applied Biosystems®, 4366596) and the miRNA-specific stem-loop primers for the miRNAs-21, -34a, -29a and -192 (Applied Biosystems™, TaqMan microRNA Assays #002438, #000426, #002112, and #000491, respectively). The reaction mix was prepared according to the manufacturer's instructions for a final reaction volume of 10 μ L with 3 μ L RNA extract. The conditions for reverse transcription were set for 30 min at 16 °C followed for 30 min at 42 °C and 5 min at 85 °C.

The real-time PCR was carried out using a fluorometric thermal cycler (Qiagen, Rotor-Gene Q 3000 or Celtic Diagnostics, Corbett research PCR Rotorgene 6000). The reaction mix contained TaqMan® Fast Advanced Master Mix 1x (Applied Biosystems™, 4444557), TaqMan microRNA Assay primer 1x (Applied Biosystems™) and 1.3 μ L cDNA in a final reaction volume of 20 μ L. The PCRs were run at 95 °C for 20 s followed by 40 cycles of 1 s at 95 °C and 20 s at 60 °C. All extracted RNAs were analyzed in technical duplicates; Ct-values of these two measurements were averaged and considered a single value.

Drug Transporter Assays

Drug-transporter interactions with P-gp and MRP2/4 were studied in ciPTEC-OAT1 using calcein-acetoxymethyl (calcein-AM, 2 μ M, Life Technologies, C1430) and 5-chloromethylfluorescein diacetate (CMFDA, 1.25 μ M, Life Technologies, C7025), respectively [16], [21]. Both calcein-AM, a P-gp substrate, and CMFDA are permeable to the cell membrane and are intracellularly metabolized into calcein or glutathione-methylfluorescein (GS-MF), substrates for MRP2/4. PSC833 (10 μ M, Tocris, Bristol, UK, 4042/1) was used as model inhibitor for P-gp, included as positive control for inhibition in each experiment. Transport of calcein and GS-MF is by MRP2/4, P-gp and BCRP and therefore, a mixture of PSC833 (10 μ M), MK571 (10 μ M, Sigma-Aldrich, M7571) and KO143 (10 μ M, Sigma-Aldrich, K2144) was used to selectively inhibit efflux, again as positive control [21]. Stocks of calcein-AM, CMFDA, PSC833 and KO143 were dissolved in DMSO (Sigma-Aldrich, D5879), MK571 was dissolved in milli-Q water. Test compounds were dissolved as described in Table S1. Work solutions were prepared in freshly-prepared Krebs-Henseleit buffer (Sigma-Aldrich, K3753), supplemented with 10 mM HEPES (Sigma-Aldrich, H3375) at pH 7.4, referred to as KHH. All conditions contained a final concentration of DMSO of 0.6% v/v. Perfusion channels were washed with 50 μ L in the medium-channel inlet and outlet and then perfused twice with KHH. Cells were incubated with calcein-AM and CMFDA with or without the presence of model inhibitors or nephrotoxicants at 37 °C for 1 h. Next, efflux was arrested using PSC833 (10 μ M), MK571 (10 μ M) and KO143 (10 μ M) in cold (4 °C) KHH. Perfusion during washing and incubation was ensured by adding 80 μ L or 20 μ L to the medium-channel inlet and outlet, respectively.

Intracellular accumulation of calcein and GS-MF in ciPTEC-OAT1 in the OrganoPlate was measured *in situ* on a spinning disk confocal Becton Dickinson (BD) Pathway 855 high-throughput microscope (BD Bioscience, Breda, The Netherlands). A $\times 10$ objective was used with excitation and emission filters set at 488 nm and 520 nm, respectively. Subsequently, bright-field images of each chip were acquired. Fluorescence intensity was determined at the ECM-medium interface in Image J. Intensity per chip was then normalized to fluorescence measured in vehicle in same experiment.

Data Analysis

Data are presented as mean \pm SD for at least three chips ($n = 3-5$) per condition, unless stated otherwise. Statistics were performed using GraphPad Prism version 8 (San Diego, CA, USA). Data was found to be significantly different compared to medium control or corresponding vehicle if $p < 0.05$ using a one-way ANOVA followed by a Tukey's multi comparison post hoc test.

Results

Robustness of the System

Cisplatin, tenofovir, cyclosporin A and tobramycin have previously been shown to cause mild to severe cytotoxicity in confluent monolayers of ciPTEC-OAT1 in a 96-well plate [14]. In our work, cell viability was measured upon a 48-h exposure to cisplatin (5 and 30 μ M), tenofovir (15.6 and 1000 μ M), cyclosporin A (5 and 30 μ M), and tobramycin (7.5 and 15 mM) in ciPTEC-OAT1 in the OrganoPlate at three different research sites. At the chosen concentrations, these nephrotoxics caused a significant reduction in cell viability (Fig. 2) and increased LDH release for tobramycin and cyclosporin A at the highest concentrations tested (Figure S2). As depicted in Fig. 2, we observed similar effects in viability of ciPTEC-OAT1 at all three laboratories, in line with earlier findings reported by lab B [14]. The other PTEC line used in the Nephroscreen, RPTEC, showed decreased viability, assessed with WST-8 and LDH-release, after exposure to tobramycin only, thus appeared to be less sensitive than ciPTEC-OAT1 (Figure S2). ciPTEC-OAT1 exposed to tenofovir in the compound screen showed a higher viability compared to the vehicle control when exposed to the low concentration but a significant loss in viability at 1000 μ M (Fig. 2).

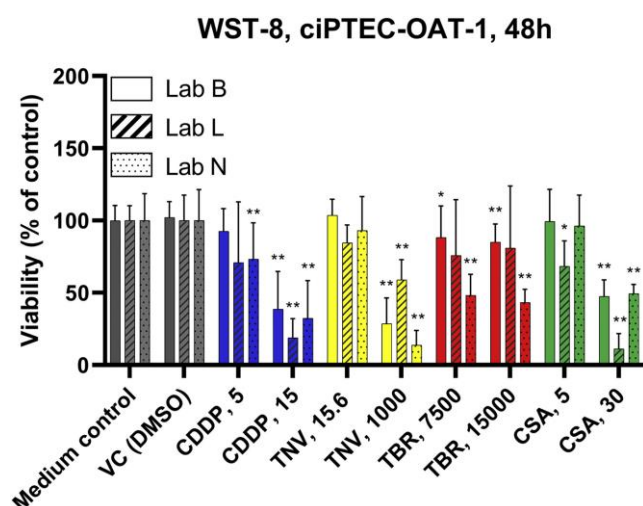


Figure 2: Side to side comparison of viability data received from the WST-8 assay. Data are from three different labs, performed at FHNW (Lab B), Mimetas (lab L), or at department of Pharmacology and Toxicology, Radboud University Medical Center (lab N). In each site, 1–3 independent experiments with 2–6 replicates (chips) were analyzed. Abbreviations: CDDP, cisplatin; CSA, cyclosporin A; DMSO, dimethyl sulfoxide; TBR, tobramycin; TNV, tenofovir; WST-8, 2-(2-methoxy-4-nitrophenyl)-3-(4-nitrophenyl)-5-(2,4-disulfophenyl)-2H-tetrazolium, monosodium salt; VC, vehicle control. #: in lab L slightly different concentrations for tenofovir were used, 15 μ M instead 15.6 μ M and 1215 μ M instead of 1000 μ M. Statistically significant compared to corresponding vehicle: * $p < 0.05$, ** $p < 0.01$.

Assessment of Cytotoxicity and Tubular Damage

Based on the data generated with the model nephrotoxics, we performed experiments to ascertain the potential nephrotoxicity of eight test substances. These compounds were provided by the pharmaceutical companies GlaxoSmithKline (GSK), Roche and Pfizer in a blinded manner,

referred to as unknown nephrotoxics, and named as G1, G2 and G3 for compounds provided by GSK, R1, R2 and R3 for compounds provided by Roche, and P1 and P2 for compounds provided by Pfizer.

In a first step, cell viability in confluent monolayers of ciPTEC-OAT1 in 2D was assessed to determine appropriate starting concentrations of the unknown nephrotoxics in the OrganoPlate (data not shown). Cytotoxicity parameters provided a consistent picture and showed reduced cell viability (WST-8) and increased LDH release caused by five of the unknown nephrotoxics (G1, G3, R1, R3, P2) in both cell lines tested (Fig. 3a–h). NAG release was only determined in ciPTEC-OAT1 and yielded concordant results (Fig. 3i and j). In general, toxicity was more marked after 48 h than after 24 h exposure. Interestingly, G2 increased cellular production of formazan in RPTECs, while R2 increased this production in ciPTEC-OAT1. The sensitivity of both PTEC cell lines to the treatments was, however, comparable.

In addition to cytotoxicity measurements, the impact of treatments on the functionality of the tubular epithelial layer (i.e. barrier integrity in RPTEC in the OrganoPlate) was determined by assessing the apparent permeability to the low-molecular weight marker dextran-TRITC (4.4 kDa) and high-molecular weight marker dextran-FITC (155 kDa). Impaired barrier function was observed after exposure to tobramycin, and five test substances (G1, G3, R1, R3 and P2). The latter were the same five compounds that caused a decrease in cell viability. The effects were generally more marked after 48 h than after 24 h exposure (Fig. 4).

Drug-Transporter Interactions

Drug-transporter interaction with P-gp and MRP2/4 is another important feature that can result in high intracellular concentrations in PTEC leading to renal toxicity of substances. Drug-transporter interactions were determined in ciPTEC-OAT1 as stable expression of these drug transporters has been previously confirmed in this model [16]. For MRP2/4, interactions were found for cyclosporin A, G2, G3 and P2 (Fig. 5a and c). Cyclosporin A, G3, P1 and P2 resulted in interactions for P-gp (Fig. 5b and d). It is interesting to note that the compound G2 clearly impaired MRP2/4 transport but did not cause any cytotoxicity based on other assays performed, suggesting competition for transport solely.







miRNA-Release

Released miRNAs can act as biomarkers of toxicity that are detectable in cell culture medium before other toxicity biomarkers. The medium of the cells treated with the eight blinded compounds was collected after 24 h of exposure and levels of four selected miRNAs determined as described in materials and methods. Most tested compounds led to increases in the miRNA panel at both tested concentrations (Fig. 6a–d). Exceptions to this assessment were compounds G2 and R2 that did not cause significantly increased release of any of the tested miRNA into the medium. Interestingly, compound R3 led to an increase in all four miRNAs at the lowest tested


concentration (25 μM) while the cells treated with 100 μM did not show changes in miRNA release.

Specific Exposures to Compounds R1 and R2


Due to the lack of toxicity observed in most parameters with compound R2 and to the high concentrations required to elicit a response with R1, tailor-made subsequent experiments were planned in discussion with the sponsors, who knew the identity of the compounds. To this end, ciPTEC-OAT1 in the OrganoPlate were exposed to several concentrations (ranging between 15 and 1000 μM) of R1 and R2 for 48 h and 11 days, respectively. At 48 h, compound R1 showed a dose dependent loss in viability assessed with the WST-8 assay with an EC50 of 214.5 and 277.3 μM (results from two independent experiments). This toxicity was accompanied by a dose dependent induction of HMOX1 (Fig. 7a and b). The results confirm that concentrations $\geq \sim 100$ μM lead to toxicity with an EC50 214 and 277 μM (two independent experiments), whereas HMOX1 expression levels increased in a concentration-dependent manner from a concentration of 62.5 μM (≥ 2 -fold induction; Fig. 7). R1 did not lead to any change in expression of NGAL or in NAG released into the medium (Fig. 7b and c). Based on information provided by the sponsors, compound R2 was additionally tested after long-term (11 days) of incubation. Cytotoxicity assay at 11 days showed a decrease in cell viability with an EC50 of 367.5 and 267.4 μM (results from two independent experiments) (Fig. 8a). As with compound R1, HMOX1 expression increased in a dose dependent manner from a concentration of 30.25 μM (≥ 2 -fold induction; Fig. 8b). The compound, however, did not cause an increase in NGAL expression or led to a release of NAG release into the medium at any of the tested time points (48 hours or 11 days of exposure, Fig. 7b and c).

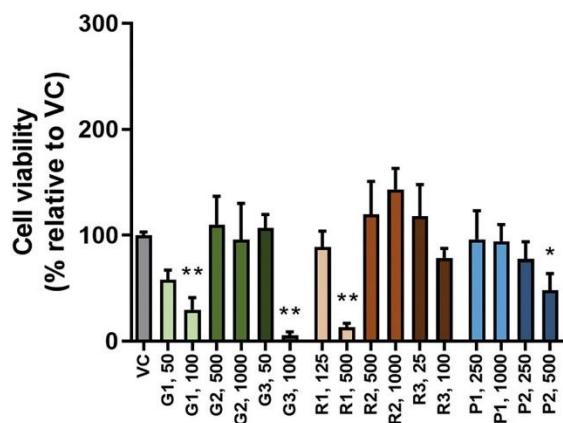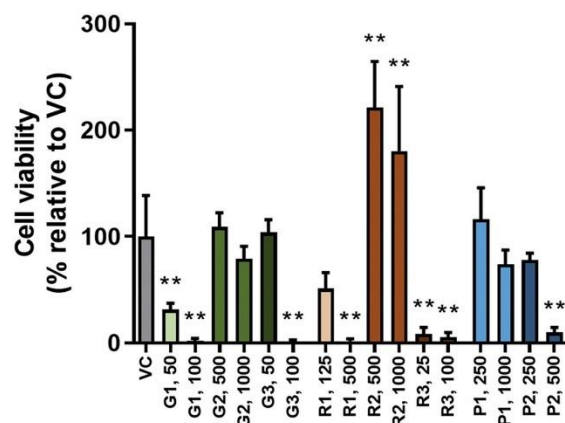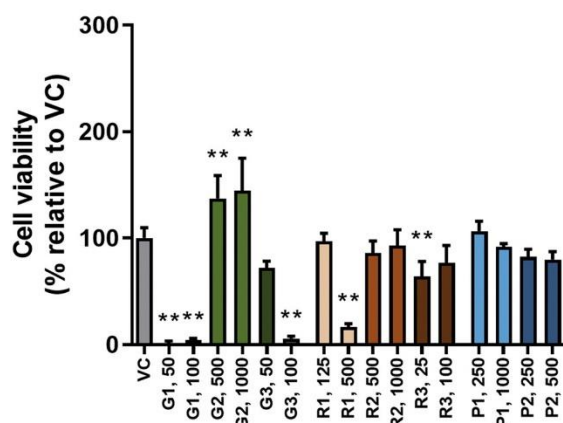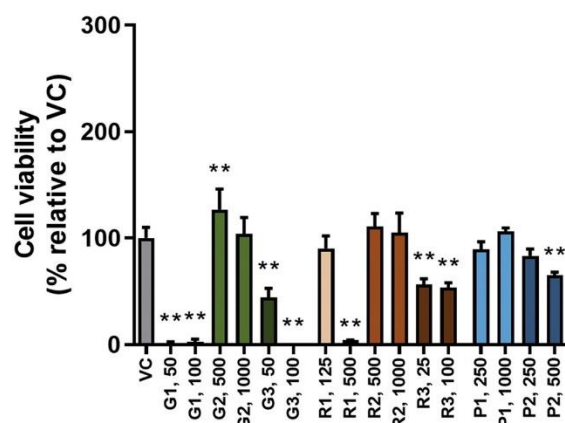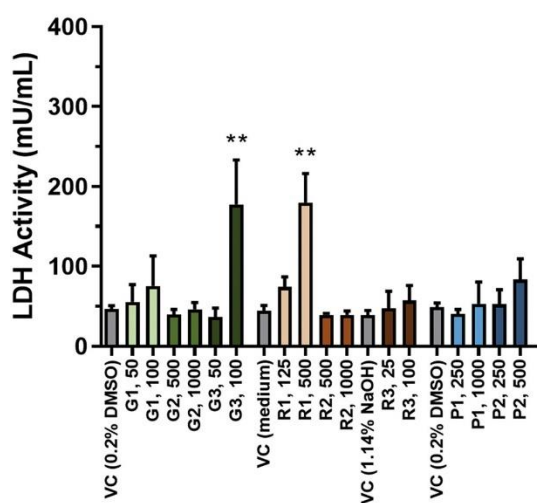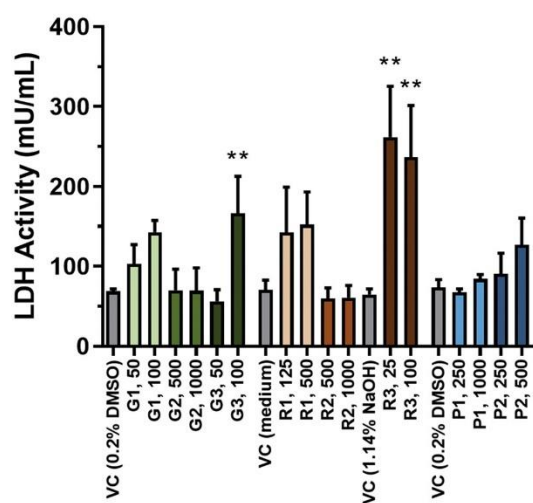a**WST-8, ciPTEC-OAT1, 24h****b****WST-8, ciPTEC-OAT1, 48h****c****WST-8, RPTEC, 24h****d****WST-8, RPTEC, 48h****e****LDH activity, ciPTEC-OAT1, 24h****f****LDH activity, ciPTEC-OAT1, 48h**


g

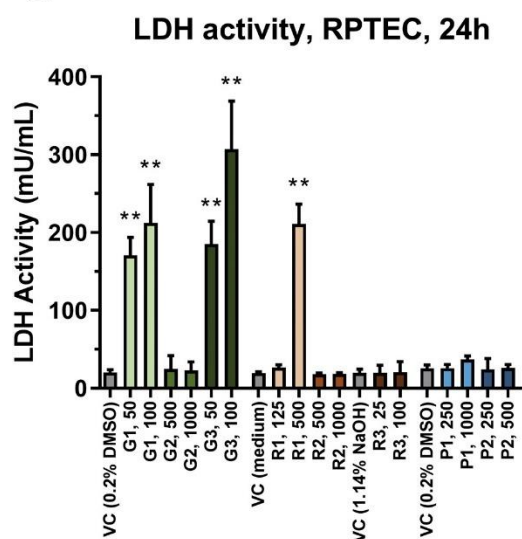

h

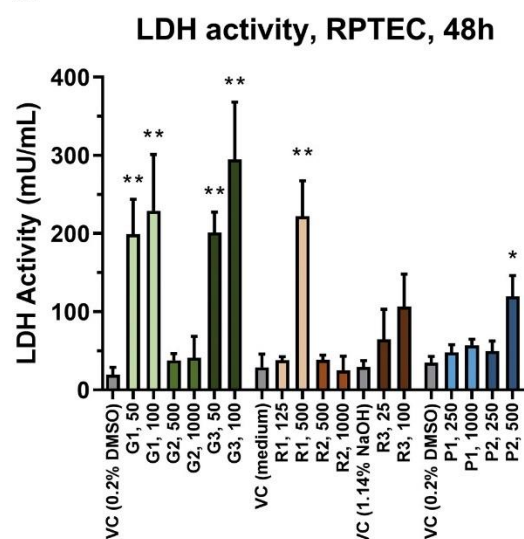

i

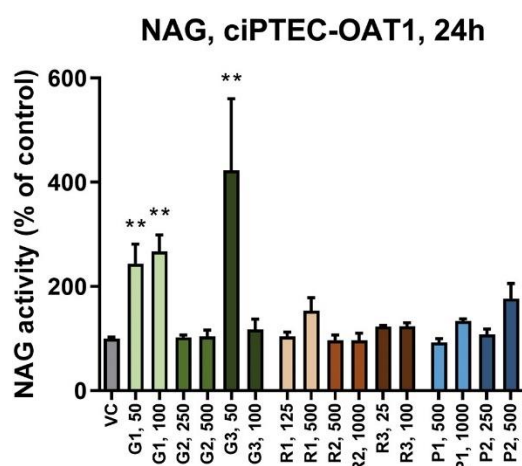

j



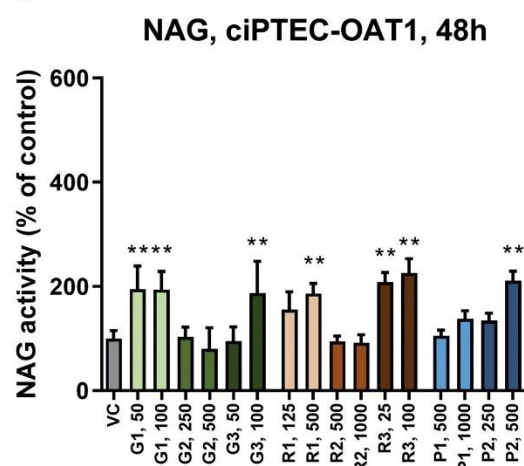


Figure 3: Cell viability and lactate dehydrogenase (LDH) release in compound screen. (a) Cell viability in ciPTEC-OAT1 after 24 h and (b) 48 h exposure. (c) Cell viability in RPTEC after 24 h and (d) 48 h exposure. (e) Release of LDH, as measure of membrane integrity, in ciPTEC-OAT1 after 24 h or (f) 48 h exposure. (g) LDH release in RPTEC after 24 h or (h) 48 h exposure. (i) N-acetyl- β -d-glucosaminidase (NAG) release by ciPTEC-OAT1 treated with sponsor compounds screen after 24 h (i) and 48 h (j) exposure. Abbreviations: VC, vehicle control; WST-8, 2-(2-methoxy-4-nitrophenyl)-3-(4-nitrophenyl)-5-(2,4-disulfophenyl)-2H-tetrazolium, monosodium salt; DMSO, Dimethyl sulfoxide; NaOH, Sodium hydroxide. Statistically significant compared to corresponding vehicle: * $p < 0.05$, ** $p < 0.01$.


a Barrier Integrity, 4.4 kDa, 24h

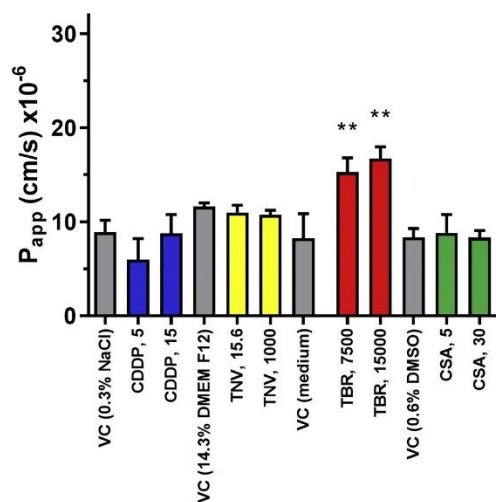


b Barrier Integrity, 4.4 kDa, 48h

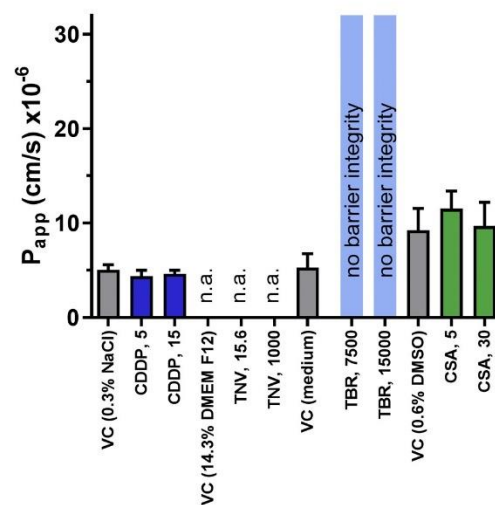


c Barrier Integrity, 150 kDa, 24h

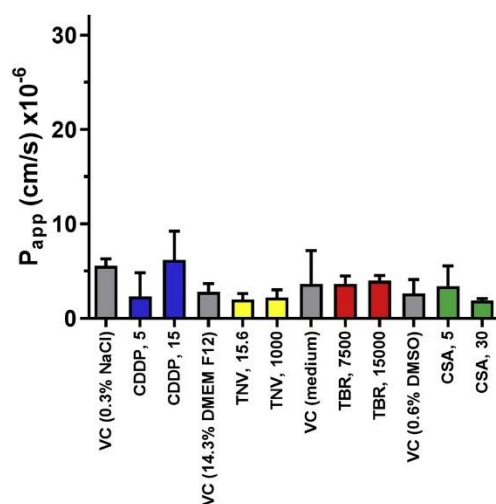


d Barrier Integrity, 150 kDa, 48h

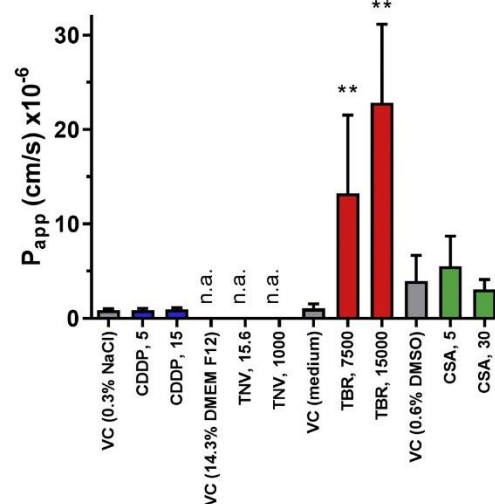

e Barrier Integrity, 4.4 kDa, 24h

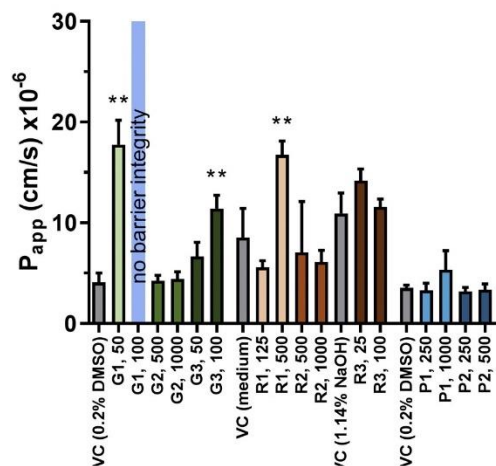


f Barrier Integrity, 4.4 kDa, 48h

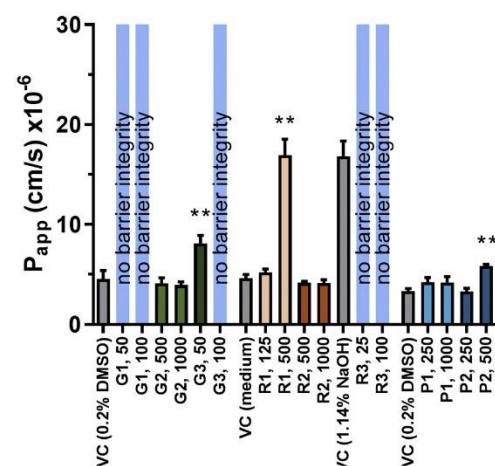

g

Barrier Integrity, 150 kDa, 24h

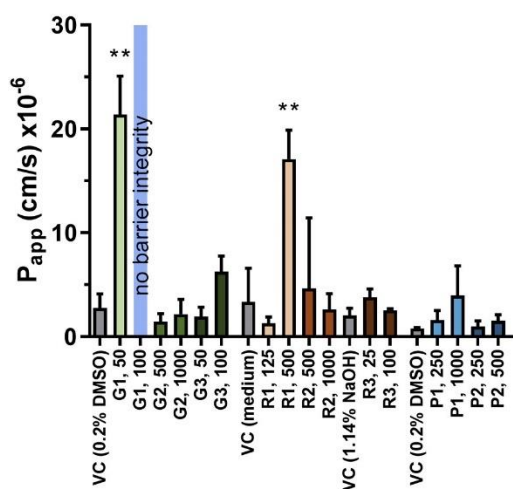

h

Barrier Integrity, 150 kDa, 48h

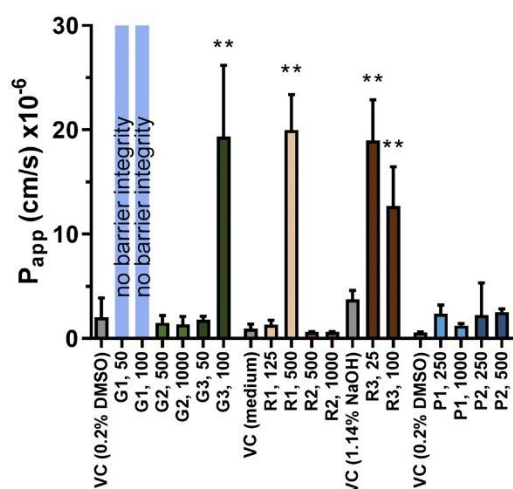


Figure 4: Barrier integrity of RPTEC tubules. (a) Leakage of high-weight molecular marker (155 kDa) fluorescein isothiocyanate (FITC)-dextran after 24 h or (b) 48 h of exposure to model nephrotoxicants. (c) Leakage of low-weight molecular marker (4.4 kDa) tetramethylrhodamine (TRITC)-dextran after 24 h or (d) 48 h of exposure to model nephrotoxicants cisplatin (CDDP), tenofovir (TNV), tobramycin (TBR) and cyclosporin A (CSA). (e) Leakage of FITC-dextran after 24 h and (f) 48 h of exposure to unknown nephrotoxicants. (g) Leakage of TRITC-dextran after 24 h or (h) 48 h of exposure to unknown nephrotoxicants. Abbreviations: D-F12, Dulbecco's modified Eagle's medium and nutrient mixture F-12; VC, vehicle control; P_{app}, apparent permeability; DMSO, Dimethyl sulfoxide; NaOH, Sodium hydroxide; kDa, kilodaltons; n.a., not available. Statistically significant compared to corresponding vehicle: *p < 0.05, **p < 0.01.

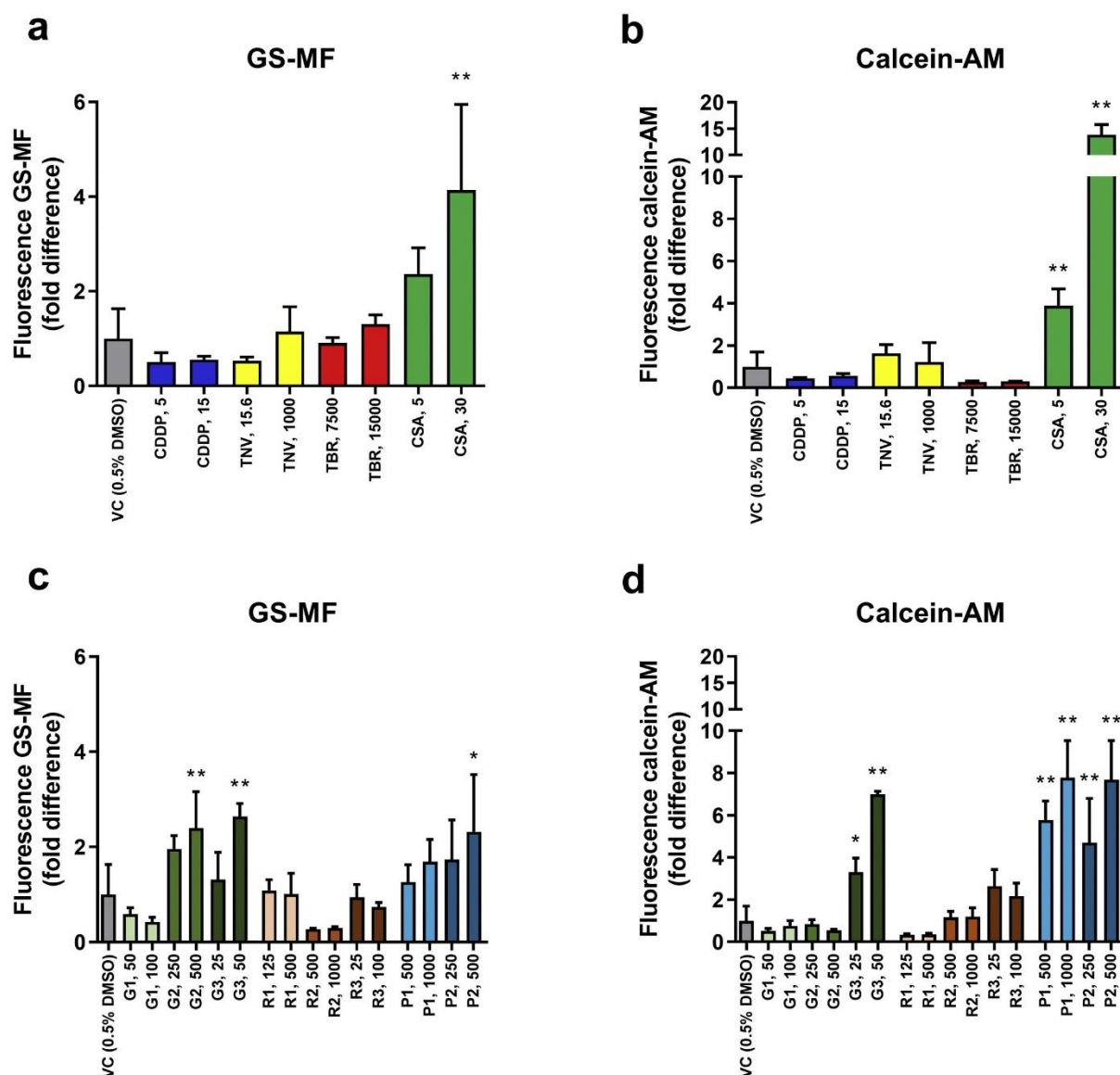


Figure 5: Drug-transporter interaction at multidrug resistance protein 2/4 (MRP2/4) and P-glycoprotein (P-gp) in ciPTEC-OAT1 after co-incubation with compounds at 37 °C for 1 h followed by arresting of efflux with PSC833, MK571 and KO143 (all 10 μ M). For model nephrotoxics cisplatin (CDDP), tenofovir (TNV), tobramycin (TBR) and cyclosporin A (CSA) (a) interactions at MRP2/4 and (b) P-gp and for unknown nephrotoxics (c) interactions at MRP2/4 and (d) P-gp. Abbreviations: GS-MF, Glutathione methylfluorescein; AM, acetoxymethyl; VC, vehicle control; DMSO, Dimethyl sulfoxide. Statistically significant compared to corresponding vehicle: * $p < 0.05$, ** $p < 0.01$.

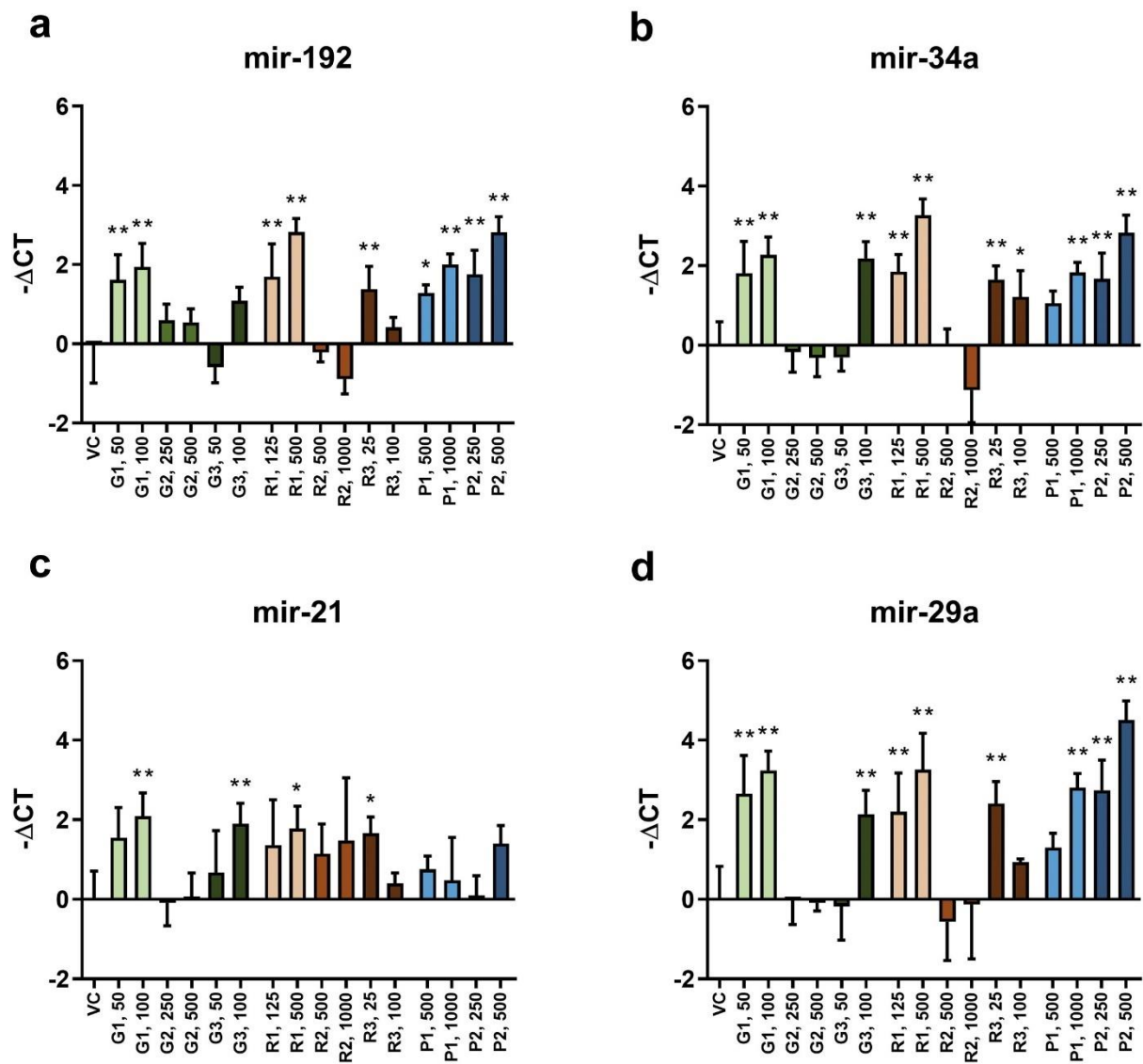


Figure 6: Release into the medium of mir-192 (a), mir-34a (b), mir-21 (c) and mir-29a (d) after 24 h exposure of ciPTEC-OAT1 to sponsor compounds. Data are represented as Delta CT values in comparison with the vehicle control (VC). Statistically significant compared to corresponding vehicle: * $p < 0.05$, ** $p < 0.01$.

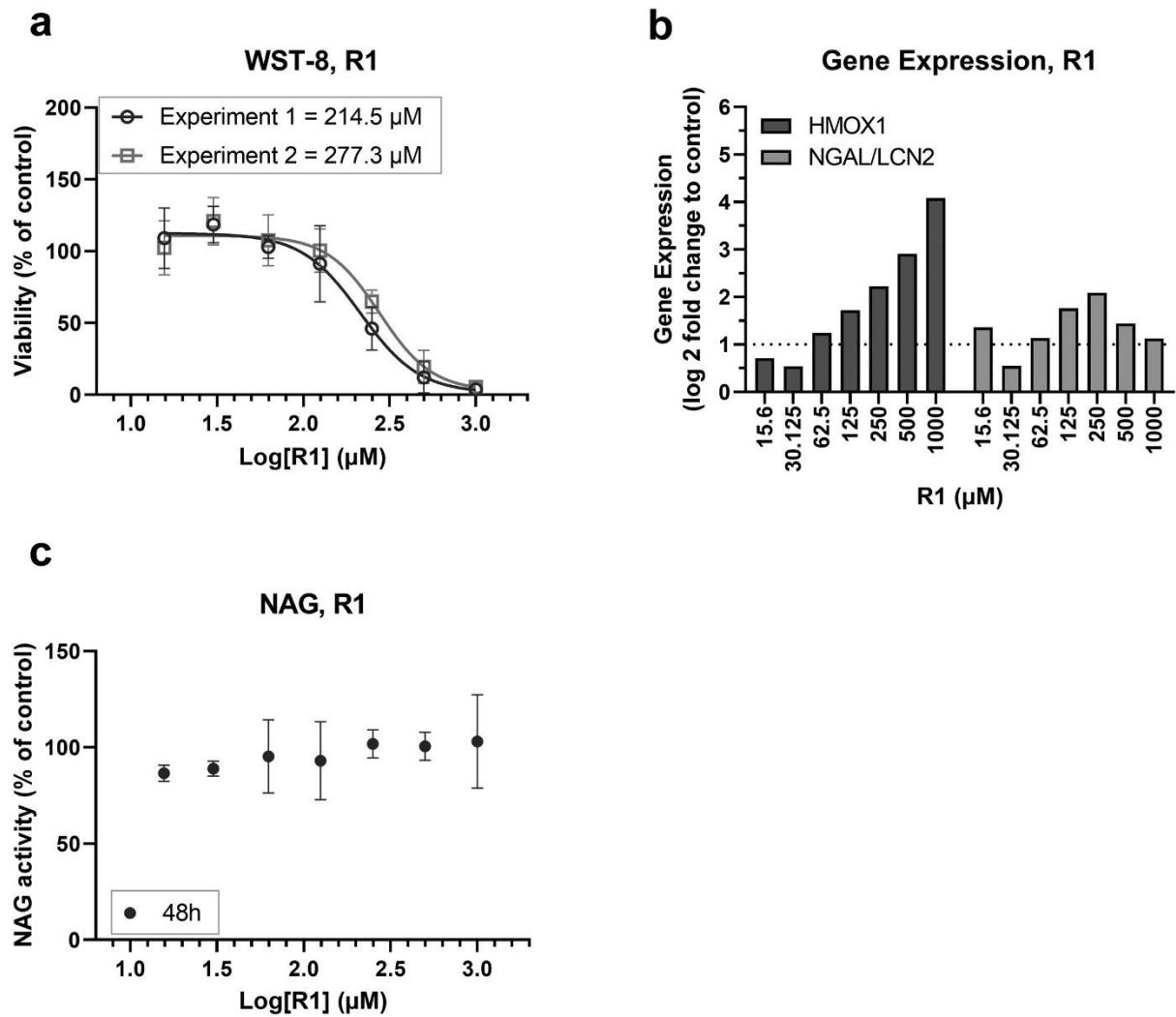


Figure 7: Exposure of *ciPTEC-OAT1* to R1 (Colistin). Data represent cell viability with the IC_{50} value of each experiment. (a), Gene expression of toxicity markers (b) and release of NAG into the medium after 48 h (c). Abbreviations: HMOX1, heme oxygenase (decycling) 1; LCN2, Lipocalin-2 (LCN2); NGAL, oncogene 24p3 or neutrophil gelatinase-associated lipocalin; NAG, N-acetyl- β -d-glucosaminidase.

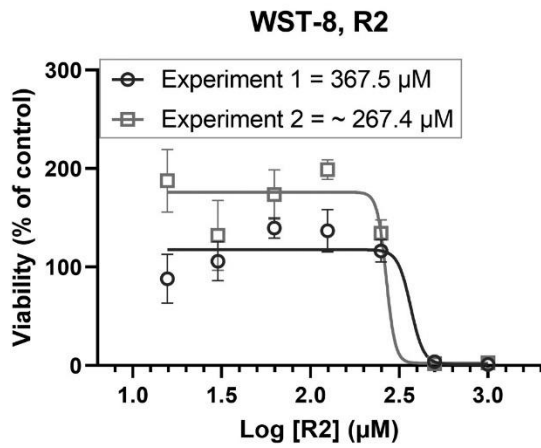
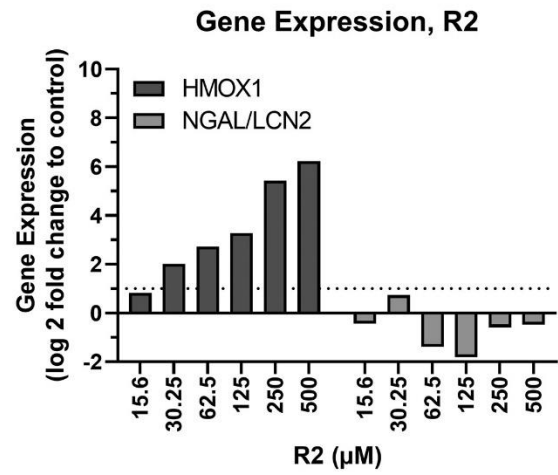
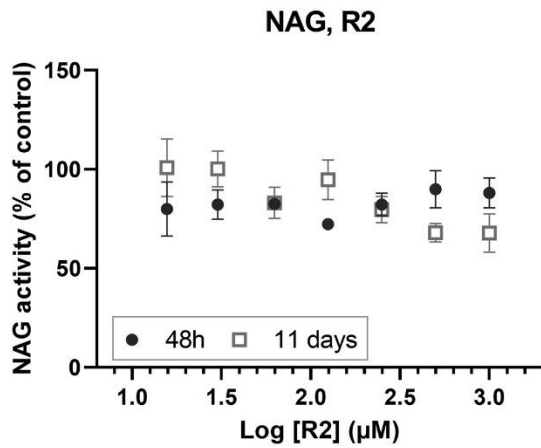
a**b****c**

Figure 8: Exposure of ciPTEC-OAT1 to R2 (Cefepime) for 48 h and 11 days. Data represent cell viability. (a), Gene expression of toxicity markers (b) and release of NAG into the medium after 48 h (circles) and 11 days (squares) (c). Abbreviations: HMOX1, heme oxygenase (decycling) 1; LCN2, Lipocalin-2 (LCN2); NGAL, oncogene 24p3 or neutrophil gelatinase-associated lipocalin; NAG, N-acetyl- β -d-glucosaminidase.

Discussion

We demonstrated the use of Nephroscreen, a proximal tubule-on-a-chip platform for the screening of nephrotoxicity and drug-transporter interactions. The robustness and transferability of this platform, a key requirement for broad implementation in industry, was established by the highly comparable results obtained at three different sites as depicted in Fig. 2. Two human renal cell lines, ciPTEC-OAT1 and RPTEC, were exposed to two concentrations of each substance (four model nephrotoxics and eight blinded compounds) at two exposure times (24 and 48 h). Subsequently, tailor-made assays were performed for two of the compounds (R1, colistin and R2, cefepime), exposing ciPTEC-OAT1. Experimental procedures were carried out in three different laboratories as depicted in Fig. 1, following established standard operating procedures (SOPs).

Concordant with the goal of the study, Nephroscreen was specifically designed to fulfill requirements of pharmaceutical companies. In pre-clinical toxicity studies, multiple drug candidates, typically up to 200 compounds, are assessed *in vitro* and/or *in vivo* in animal experimentation [22]. In this study, we implemented well-characterized PTEC lines, considered relevant for nephrotoxicity and drug-transporter interaction studies, as potential alternatives to animal experimentation [14], [15], [16], [17], [18], [23], [24]. Supply and reproducibility (low batch-to-batch variability) of these commercially available cells are guaranteed. The PTECs were cultured in the OrganoPlate, generating a proximal tubule-on-a-chip consisting of 40 chips on a 384-well microtiter plate format [15]. This model enables the culture of the cells as a tubular structure kept under flow, without the use of support membranes. Perfusion was generated by passive leveling resulting in a bidirectional, oscillating flow allowing significant levels of shear stress circumventing the use of pumps [15], [25]. The diameter of the tubule (approximately 400 μm) is larger than in the human proximal tubule. However, it represents a useful model that allows the long-term culture of polarized cell layers against the ECM with significant barrier function [15].

Effect of the treatments on the PTECs was determined by measuring functional and biochemical parameters. Functional parameters focused on the capacity of compounds to interact with transporters and to impair barrier function. Membrane drug transporters are an important characteristic of PTECs and expression and activity of transporters have been demonstrated previously [6], [17], [18], [24]. Differences in response towards toxicant exposures between ciPTEC-OAT1 and RPTEC could be explained by different expression levels of drug transporters. For instance, ciPTEC-OAT1, but not RPTEC, functionally express OAT1. This explains the lack of sensitivity of RPTEC towards tenofovir, as this substance requires the OAT1 transporter for cellular uptake [18], [24]. Epithelial barrier function is also a key functional parameter, as a leaky epithelium is often an indicator of impaired kidney function, previously demonstrated in this model [15].

Regarding biochemical parameters, cell viability, LDH release and gene expression of toxicity markers were determined on both PTEC lines to assess the effect of the compounds. For ciPTEC-OAT1, assays for drug-transporter interaction, NAG release and miRNA secretion were also performed [14],[16], while RPTEC tubules were more suitable for the assessment of barrier integrity [15]. Cell viability, as measured by enzymatic production of formazan (WST-8 assay), and cell membrane integrity, determined by LDH release into the medium, performed well as *in vitro* markers for cytotoxicity, although with different kinetic profiles (Table 1, Figure S2). In ciPTECs, leakage of NAG into the medium displayed a slightly higher sensitivity in detecting the toxicity caused by compounds G1 and P2 than LDH release (Fig. 3), but a larger data set would be necessary to corroborate this claim. Decreased barrier integrity was almost always coupled with decreased cell viability and increased LDH release. Thus, we assume that the loss of barrier integrity was directly associated with increased cytotoxicity (Table 1).

The four model nephrotoxicants (cisplatin, tenofovir, tobramycin, cyclosporin A) assessed in this study elicited a response in at least four of the end-point read-outs (Figure S2, Table 1), in agreement with previous results [26], [27], [28], [29], [30]. Furthermore, our data show that cyclosporin A interacts with P-gp, corroborating that Nephroscreen results reflect its known P-gp inhibition [31]. Nephrotoxicity of tobramycin, on the other hand, could only be observed at concentrations that were up to 1500-fold higher than plasma concentrations generally found in patients (0.01 mM) [30]. This is probably due to the low expression of cubilin and megalin receptors in ciPTEC-OAT1 [17], [32], [33] a limitation that needs to be taken into consideration for compounds known to be dependent on this cellular uptake mechanism. During the second phase of the study, eight substances selected and provided by the sponsors were assessed in a blinded manner.

1000 μ M and unpublished *in vitro* results with PTEC toxicity at concentrations >500 μ M (Table S2).

The ideal set-up for a screening tool requires short-term incubation in order to generate results quickly. Our results show that for most compounds, exposure during 24 and 48 h suffices to detect compound-induced damage to renal tubular cells. However, as demonstrated with the results obtained with R2, the developed platform can also be implemented as a second tier assay for selected substances of interest that may require subchronic or chronic exposure. The sensitivity and specificity of Nephroscreen cannot be determined as only a small set of test compounds was included. Future testing should expand the number of compounds and include compounds with other toxicity target organs (non-nephrotoxic) as well as non-toxic compounds to further characterize the model. Also, side-by-side comparison with conventional 2D cell cultures may be performed to achieve direct comparison and therefore show the value in this model over more simple models. Other proximal tubule-on-a-chip models have, however, demonstrated that 3D microfluidic models increase *in vivo* physiology and sensitivity of PTECs towards nephrotoxicant

exposure [11], [12], [13]. Additional future improvements could include the corroboration and systematic assessment of HMOX-1 induction as an earlier biomarker and its implementation as a FRET-based assay for HMOX1 expression [19]. The implementation of cell lines expressing more types of functional transporters relevant, such as ciPTEC-OAT1/OAT-3 could also be envisaged [18]. Moreover, data on additional compounds could support the creation of intelligent decision algorithms to evaluate the multiparametric data provided by Nephroscreen in an unbiased manner.

Preliminary knowledge such as clinical data (for R1, R2, R3) or *in vivo* toxicity data were provided by sponsors after the Nephroscreen data collection and analysis had been finalized. The compound concentrations selected and applied to Nephroscreen unveiled nephrotoxic potential at concentrations that are considered relevant for toxicological assessment (Table S2). For R1, R2, and R3 there is clinical data available; these compounds are used at very high concentrations and administered parenterally, so that a high systemic exposure is reached in patients. For the other blinded compounds, there is no clinical data but preclinical data on at least one animal species show that the exposure (C_{max}) at which nephrotoxicity was observed in a subset of the animals was slightly lower, but in a similar range to that eliciting a positive flag in Nephroscreen. In our experimental set-up, prediction of nephrotoxic liability was highest when combining results of cell viability, LDH release, and miRNA release. Additional functional effects of the compounds were uncovered by the drug-transport interaction measurements. The majority of the unknown nephrotoxicants (G1, G3, R1, R2, R3, P1, and P2) resulted in toxicity detected by at least one of the parameters or in drug-transporter interactions. Interestingly, an interaction at MRP2/4 was observed for G2, despite the lack of toxicity. This could be explained by the fact that G2-induced nephrotoxicity was only observed upon long-term exposure *in vivo* (one month, highest tested dose) as compared to short-term (up to 48 h) exposure tested in Nephroscreen. The compound P1 showed a particularly mild effect on Nephroscreen, as it only led to increased levels of miRNAs (mir192, mir29a, and mir34a) in the medium and an interaction at P-gp, without causing cytotoxicity. This is consistent with the existing animal data showing that this vasopressin 1-a antagonist caused tubular degeneration/regeneration only in female rats treated with a high dose (125 mg/kg/day) for two weeks.

Table 1: Overview Read-Outs Nephrotoxicity and Renal Drug-Transporter Interactions Assessed in This Nephrotoxicant Screening.

	Conc. (µM)	ciPTEC-OAT1			Drug-Transporter Interaction	NAG Activity	miRNA Release	RPTEC		
		Cell Viability	LDH Release					Cell Viability	LDH Release	Barrier Integrity
CDDP	5									
	15	↓ 48 h								
TNV	15.6	↑ 48 h								
	1000	↓ 48 h								
TBR	7500									
	15000		↑ 24 h					↓ 48 h	↑ 48 h	4.4 kDa: ↓ 24 h 150 kDa: ↓ 48 h
CSA	5				P-gp					
	30	↓ 48 h	↑ 24 h		MRP2/4, P-gp					
G1	50	↓ 48 h				↑ 24 h	mir-34a, mir-29a, mir-192	↓ 24 h	↑ 24 h	4.4 kDa: ↓ 24 h 150 kDa: ↓ 24 h
	100	↓ 24 h				↑ 24 h	mir-21, mir-34a, mir-29a, mir-192	↓ 24 h	↑ 24 h	4.4 kDa: ↓ 24 h 150 kDa: ↓ 24 h
G2	250									
	500				MRP2/4			↑ 24 h		
	1000							↑ 24 h		
G3	25				P-gp					
	50				MRP2/4, P-gp	↑ 24 h	mir-21, mir-34a, mir-29a	↓ 48 h	↑ 24 h	4.4 kDa: ↓ 48 h
	100	↓ 24 h	↑ 24 h			↑ 48 h	mir-34a, mir-29a, mir-192	↓ 24 h	↑ 24 h	4.4 kDa: ↓ 24 h 150 kDa: ↓ 48 h
R1	125									
	500	↓ 24 h				↑ 48 h	mir-21, mir-34a, mir-29a, mir-192	↓ 24 h	↑ 24 h	4.4 kDa: ↓ 24 h 150 kDa: ↓ 24 h
R2	500	↑ 48 h								
	1000	↑ 48 h								
R3	25	↓ 48 h	↑ 48 h			↑ 48 h	mir-21, mir-34a, mir-29a, mir-192	↓ 24 h		4.4 kDa: ↓ 48 h 150 kDa: ↓ 48 h
	100	↓ 48 h	↑ 48 h			↑ 48 h	mir-34a	↓ 48 h		4.4 kDa: ↓ 48 h 150 kDa: ↓ 48 h
P1	250						mir-192			
	500				P-gp					
	1000				P-gp		mir-34a, mir-29a, mir-192			
P2	250				P-gp		mir-34a, mir-29a, mir-192			
	500	↓ 24 h			MRP2/4, P-gp	↑ 48 h	mir-34a, mir-29a, mir-192	↓ 48 h	↑ 48 h	4.4 kDa: ↓ 48 h

For cell viability, LDH release, drug-transporter interaction, and barrier integrity: ↑ represents significant increased effect observed compared to corresponding vehicle, ↓ represents significant decreased effect observed compared to corresponding vehicle. CDDP, cisplatin; CLDN2, claudin-2; Conc., concentration; CSA, cyclosporin A; FITC, fluorescein isothiocyanate; LDH, lactate dehydrogenase; MRP2/4, multidrug resistance-associated proteins 2 and 4; P-gp, P-glycoprotein; TBR, tobramycin; TNV, tenofovir; TRITC, tetramethylrhodamine.

In the last phase of this study, tailor-made experiments were performed including expanded concentration ranges of two compounds after their identities had been uncovered: R1 (colistin) and R2 (cefepime). For R1 (colistin) we repeated the cytotoxicity assay and determined gene expression of HMOX-1 and NGAL, as well as the release of NAG in the medium. The EC₅₀ for cytotoxicity obtained in this experiment was in line with the results of the first round of experiments (~200 μ M). The *in vitro* cytotoxicity of colistin has been reported by others in the high micromolar range, around 690 μ M.³⁴ Preliminary gene expression data, however, uncovered a dose-dependent increase in HMOX-1 from a concentration of 30 μ M (>2-fold induction). Transcriptional upregulation of heme oxygenase-1 is a well-known indicator of cellular stress in several cells and tissues, including PTECs and may be a more sensitive biomarker of toxicity [19], [26]. This was also supported by the long-term (11 days) exposure of the Nephroscreen to R2 (cefepime). Cefepime had an EC₅₀-value for cytotoxicity of approximately 300 μ M, while the concentration-dependent transcriptional induction of HMOX-1 was observed at concentrations from 62 μ M (>2-fold induction). These results also show that a long-term exposure is required for this compound, as the short-term experiment failed to show toxicity up to a concentration of 1000 μ M and unpublished *in vitro* results with PTEC toxicity at concentrations >500 μ M (Table S2).

The ideal set-up for a screening tool requires short-term incubation in order to generate results quickly. Our results show that for most compounds, exposure during 24 and 48 h suffices to detect compound-induced damage to renal tubular cells. However, as demonstrated with the results obtained with R2, the developed platform can also be implemented as a second tier assay for selected substances of interest that may require subchronic or chronic exposure. The sensitivity and specificity of Nephroscreen cannot be determined as only a small set of test compounds was included. Future testing should expand the number of compounds and include compounds with other toxicity target organs (non-nephrotoxic) as well as non-toxic compounds to further characterize the model. Also, side-by-side comparison with conventional 2D cell cultures may be performed to achieve direct comparison and therefore show the value in this model over more simple models. Other proximal tubule-on-a-chip models have, however, demonstrated that 3D microfluidic models increase *in vivo* physiology and sensitivity of PTECs towards nephrotoxicant exposure [11], [12], [13]. Additional future improvements could include the corroboration and systematic assessment of HMOX-1 induction as an earlier biomarker and its implementation as a FRET-based assay for HMOX1 expression [19]. The implementation of cell lines expressing more types of functional transporters relevant, such as ciPTEC-OAT1/OAT-3 could also be envisaged [18]. Moreover, data on additional compounds could support the creation of intelligent decision algorithms to evaluate the multiparametric data provided by Nephroscreen in an unbiased manner.

Conclusion

An advanced screening tool is needed to increase speed, reduce costs and animal testing while assessing the potential nephrotoxicity of new drug candidates. Taken together, our results indicate that Nephroscreen, consisting of a microfluidic organ-on-a-chip system coupled with a multiparametric biomarker analysis is able to identify potential nephrotoxicants. Nephroscreen is a reliable medium-throughput, standardized, automatable system that proved efficacious in identifying nephrotoxicants and provided insights into their mode of toxicity.

Conflicts of Interest

This publication contains original work. The authors Marianne K. Vormann, Henriëtte L. Lanz, Linda Gijzen, Angelique van den Heuvel, Simon Hutter, Jos Joore, Sebastiaan J. Trietsch, and Paul Vulto are or were employees of Mimetas B.V. This affiliation is declared. The OrganoPlate® is a registered trademark of MIMETAS BV. Martijn Wilmer and Frans Russel are co-inventors on patent EP2010/066792 “Novel conditionally immortalized human proximal tubule cell line expressing functional influx and efflux transporters” assigned to Radboud University Medical Center and commercialized through Cell4Pharma.com, for which Martijn Wilmer acts as scientific advisor.

Acknowledgements

Funding: This project was supported under the Crack-it challenge 15 (Nephrotube) project no. 37497-25920, an initiative of the National Center for the Replacement, Refinement and Reduction of Animals in Research (NC3Rs), UK.

References

- [1] Soo JY, Jansen J, Masereeuw R, Little MH 2018. Advances in predictive *in vitro* models of drug-induced nephrotoxicity. *Nat Rev Nephrol* 14(6):378-393.
- [2] Uchino S, Kellum JA, Bellomo R, Doig GS, Morimatsu H, Morgera S, Schetz M, Tan I, Bouman C, Macedo E, Gibney N, Tolwani A, Ronco C, Beginning, Ending Supportive Therapy for the Kidney I 2005. Acute renal failure in critically ill patients: a multinational, multicenter study. *JAMA* 294(7):813-818.
- [3] Hoste EA, Bagshaw SM, Bellomo R, Cely CM, Colman R, Cruz DN, Edipidis K, Forni LG, Gomersall CD, Govil D, Honore PM, Joannes-Boyau O, Joannidis M, Korhonen AM, Lavrentieva A, Mehta RL, Palevsky P, Roessler E, Ronco C, Uchino S, Vazquez JA, Vidal Andrade E, Webb S, Kellum JA 2015. Epidemiology of acute kidney injury in critically ill patients: the multinational AKI-EPI study. *Intensive Care Med* 41(8):1411-1423.
- [4] Mehta RL, Pascual MT, Soroko S, Savage BR, Himmelfarb J, Ikizler TA, Paganini EP, Chertow GM, Program to Improve Care in Acute Renal D 2004. Spectrum of acute renal failure in the intensive care unit: the PICARD experience. *Kidney Int* 66(4):1613-1621.
- [5] Guengerich FP 2011. Mechanisms of drug toxicity and relevance to pharmaceutical development. *Drug Metab Pharmacokinet* 26(1):3-14.
- [6] Morrissey KM, Stocker SL, Wittwer MB, Xu L, Giacomini KM 2013. Renal transporters in drug development. *Annu Rev Pharmacol Toxicol* 53:503-529.
- [7] Varma MV, Feng B, Obach RS, Troutman MD, Chupka J, Miller HR, El-Kattan A 2009. Physicochemical determinants of human renal clearance. *J Med Chem* 52(15):4844-4852.
- [8] Ivanyuk A, Livio F, Biollaz J, Buclin T 2017. Renal Drug Transporters and Drug Interactions. *Clin Pharmacokinet* 56(8):825-892.
- [9] Wilmer MJ, Ng CP, Lanz HL, Vulto P, Suter-Dick L, Masereeuw R 2016. Kidney-on-a-Chip Technology for Drug-Induced Nephrotoxicity Screening. *Trends Biotechnol* 34(2):156-170.
- [10] Tiong HY, Huang P, Xiong S, Li Y, Vathsala A, Zink D 2014. Drug-induced nephrotoxicity: clinical impact and preclinical *in vitro* models. *Mol Pharm* 11(7):1933-1948.
- [11] Homan KA, Kolesky DB, Skylar-Scott MA, Herrmann J, Obuobi H, Moisan A, Lewis JA 2016. Bioprinting of 3D Convulated Renal Proximal Tubules on Perfusable Chips. *Sci Rep* 6:34845.
- [12] Jang KJ, Mehr AP, Hamilton GA, McPartlin LA, Chung S, Suh KY, Ingber DE 2013. Human kidney proximal tubule-on-a-chip for drug transport and nephrotoxicity assessment. *Integrative biology : quantitative biosciences from nano to macro* 5(9):1119-1129.
- [13] Sakolish C, Weber EJ, Kelly EJ, Himmelfarb J, Mouneimne R, Grimm FA, House JS, Wade T, Han A, Chiu WA, Rusyn I 2018. Technology Transfer of the Microphysiological Systems: A Case Study of the Human Proximal Tubule Tissue Chip. *Sci Rep* 8(1):14882.
- [14] Suter-Dick L, Mauch L, Ramp D, Caj M, Vormann MK, Hutter S, Lanz HL, Vriend J, Masereeuw R, Wilmer MJ 2018. Combining Extracellular miRNA Determination with Microfluidic 3D Cell Cultures for the Assessment of Nephrotoxicity: a Proof of Concept Study. *AAPS J* 20(5):86.

- [15] Vormann MK, Gijzen L, Hutter S, Boot L, Nicolas A, van den Heuvel A, Vriend J, Ng CP, Nieskens TTG, van Duinen V, de Wagenaar B, Masereeuw R, Suter-Dick L, Trietsch SJ, Wilmer M, Joore J, Vulto P, Lanz HL 2018. Nephrotoxicity and Kidney Transport Assessment on 3D Perfused Proximal Tubules. *AAPS J* 20(5):90.
- [16] Vriend J, Nieskens TTG, Vormann MK, van den Berge BT, van den Heuvel A, Russel FGM, Suter-Dick L, Lanz HL, Vulto P, Masereeuw R, Wilmer MJ 2018. Screening of Drug-Transporter Interactions in a 3D Microfluidic Renal Proximal Tubule on a Chip. *AAPS J* 20(5):87.
- [17] Adler M, Ramm S, Hafner M, Muhlich JL, Gottwald EM, Weber E, Jaklic A, Ajay AK, Svoboda D, Auerbach S, Kelly EJ, Himmelfarb J, Vaidya VS 2016. A Quantitative Approach to Screen for Nephrotoxic Compounds *In vitro*. *J Am Soc Nephrol* 27(4):1015-1028.
- [18] Schindelin J, Arganda-Carreras I, Frise E, Kaynig V, Longair M, Pietzsch T, Preibisch S, Rueden C, Saalfeld S, Schmid B, Tinevez JY, White DJ, Hartenstein V, Eliceiri K, Tomancak P, Cardona A 2012. Fiji: an open-source platform for biological-image analysis. *Nat Methods* 9(7):676-682.
- [19] Caetano-Pinto P, Janssen MJ, Gijzen L, Verscheijden L, Wilmer MJ, Masereeuw R 2016. Fluorescence-Based Transport Assays Revisited in a Human Renal Proximal Tubule Cell Line. *Mol Pharm* 13(3):933-944.
- [20] Sakolish CM, Esch MB, Hickman JJ, Shuler ML, Mahler GJ 2016. Modeling Barrier Tissues *In vitro*: Methods, Achievements, and Challenges. *EBioMedicine* 5:30-39.
- [21] Fedecostante M, Westphal KGC, Buono MF, Sanchez Romero N, Wilmer MJ, Kerkering J, Baptista PM, Hoenderop JG, Masereeuw R 2018. Recellularized Native Kidney Scaffolds as a Novel Tool in Nephrotoxicity Screening. *Drug Metab Dispos* 46(9):1338-1350.
- [22] Li S, Zhao J, Huang R, Steiner T, Bournier M, Mitchell M, Thompson DC, Zhao B, Xia M 2017. Development and Application of Human Renal Proximal Tubule Epithelial Cells for Assessment of Compound Toxicity. *Curr Chem Genom Transl Med* 11:19-30.
- [23] Nieskens TT, Peters JG, Schreurs MJ, Smits N, Woestenenk R, Jansen K, van der Made TK, Roring M, Hilgendorf C, Wilmer MJ, Masereeuw R 2016. A Human Renal Proximal Tubule Cell Line with Stable Organic Anion Transporter 1 and 3 Expression Predictive for Antiviral-Induced Toxicity. *AAPS J* 18(2):465-475.
- [24] Wilmer MJ, Saleem MA, Masereeuw R, Ni L, van der Velden TJ, Russel FG, Mathieson PW, Monnens LA, van den Heuvel LP, Levtchenko EN 2010. Novel conditionally immortalized human proximal tubule cell line expressing functional influx and efflux transporters. *Cell and tissue research* 339(2):449-457.
- [25] van Duinen V, van den Heuvel A, Trietsch SJ, Lanz HL, van Gils JM, van Zonneveld AJ, Vulto P, Hankemeier T 2017. 96 perfusable blood vessels to study vascular permeability *in vitro*. *Sci Rep* 7(1):18071.
- [26] Ozkok A, Edelstein CL 2014. Pathophysiology of cisplatin-induced acute kidney injury. *Biomed Res Int* 2014:967826.
- [27] Jafari A, Khalili H, Dashti-Khavidaki S 2014. Tenofovir-induced nephrotoxicity: incidence, mechanism, risk factors, prognosis and proposed agents for prevention. *Eur J Clin Pharmacol* 70(9):1029-1040.

- [28] Ramamoorthy H, Abraham P, Isaac B 2014. Mitochondrial dysfunction and electron transport chain complex defect in a rat model of tenofovir disoproxil fumarate nephrotoxicity. *J Biochem Mol Toxicol* 28(6):246-255.
- [29] Mingeot-Leclercq MP, Tulkens PM 1999. Aminoglycosides: nephrotoxicity. *Antimicrob Agents Chemother* 43(5):1003-1012.
- [30] Schentag JJ, Lasezkay G, Cumbo TJ, Plaut ME, Jusko WJ 1978. Accumulation pharmacokinetics of tobramycin. *Antimicrob Agents Chemother* 13(4):649-656.
- [31] Polli JW, Wring SA, Humphreys JE, Huang L, Morgan JB, Webster LO, Serabjit-Singh CS 2001. Rational use of *in vitro* P-glycoprotein assays in drug discovery. *J Pharmacol Exp Ther* 299(2):620-628.
- [32] Vriend J, Peters JGP, Nieskens TTG, Škovroňová R, Blaimschein N, Schmidts M, Roepman R, Schirris TJJ, Russel FGM, Masereeuw R, Wilmer MJ 2020. Flow stimulates drug transport in a human kidney proximal tubule-on-a-chip independent of primary cilia. *Biochim Biophys Acta Gen Subj* 1864(1):129433.
- [33] Sun Y, Goes Martini A, Janssen MJ, Garrelds IM, Masereeuw R, Lu X, Danser AHJ 2020. Megalin: A Novel Endocytic Receptor for Prorenin and Renin. *Hypertension* 75(5):1242-1250.
- [34] Roberts KD, Azad MA, Wang J, Horne AS, Thompson PE, Nation RL, Velkov T, Li J 2015. Antimicrobial Activity and Toxicity of the Major Lipopeptide Components of Polymyxin B and Colistin: Last-line Antibiotics against Multidrug-Resistant Gram-negative Bacteria. *ACS Infect Dis* 1(11):568-575.

Supplementary Information

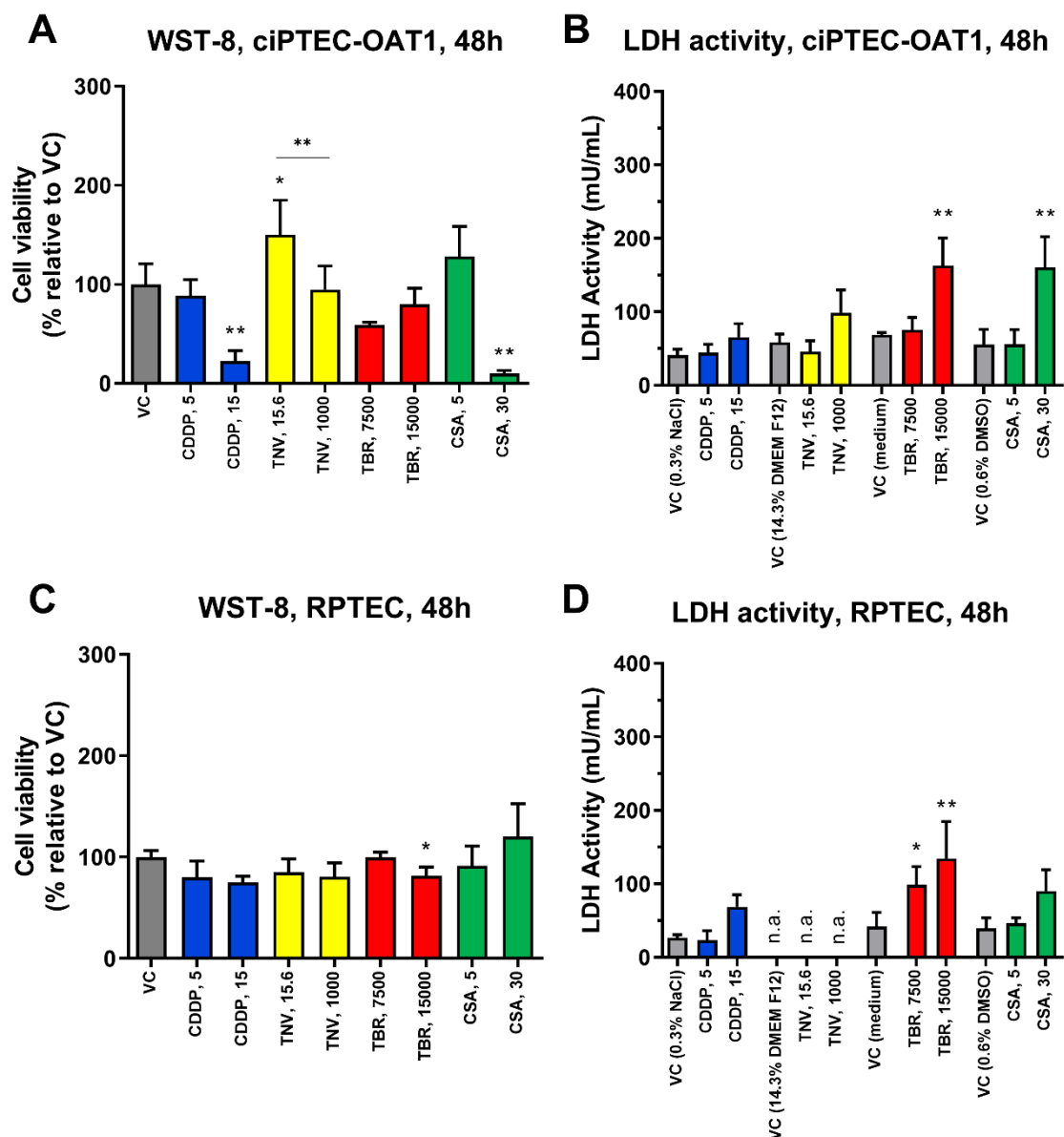


Figure S1: Release of lactate dehydrogenase (LDH) (B, D)) and Cell Viability (A, C) upon exposure to model nephrotoxics cisplatin (CDDP), tenofovir (TNV), tobramycin (TBR) and cyclosporin A (CSA) in ciPTEC-OAT1 after 48 hours. Abbreviations: D-F12, Dulbecco's modified Eagle's medium and nutrient mixture F-12; WST-8, 2-(2-methoxy-4-nitrophenyl)-3-(4-nitrophenyl)-5-(2,4-disulfophenyl)-2H-tetrazolium, monosodium salt; HMOX1, heme oxygenase (decycling) 1; LDH, lactate dehydrogenase; DMSO, Dimethyl sulfoxide; VC, vehicle control; n.a., not available. Statistically significant compared to corresponding vehicle: * $p < 0.05$, ** $p < 0.01$.

Table S1: Test compounds used to study nephrotoxicity and renal drug-transporter interactions

Compound exposure toxicity assessment										3D fluorescent efflux assays		
Compound	Supplier, catalog no.	Vehicle				Exposure				Vehicle		
		Stock conc. (mM)	Final conc. (v/v)	vehicle conc. (mM)	Exposure conc. (mM)	Stock conc. (mM)	Final conc. (v/v)	vehicle conc. (mM)	Exposure conc. (mM)	Stock conc. (mM)	Vehicle	Exposure conc. (mM)
CDDP	Sigma-Aldrich, P4394	0.9% NaCl in H ₂ O (lab L), PTECCM (lab N)	5 (lab L), 1.5 (lab N)	0.3% (lab L)	0.005, 0.015	DMSO	15	0.005, 0.015	0.005, 0.015	0.005, 0.015	0.005, 0.015	0.005, 0.015
TNV	Santa Cruz, 204335	Culture medium (lab L), HBSS (lab N)	7 (lab L), 5 (lab N)	-	0.0156, 1	HBSS	5	0.0156, 1	0.0156, 1	0.0156, 1	0.0156, 1	0.0156, 1
TBR	Sigma-Aldrich, T1783	Culture medium (lab L), PTEC SFM (lab N)	100	-	7.5, 15	KHH	100	7.5, 15	7.5, 15	7.5, 15	7.5, 15	7.5, 15
CSA	Sigma-Aldrich, 30024	DMSO	5 (lab L), 40 (lab N)	0.6% (lab L), 0.08% (lab N)	0.005, 0.03	DMSO	40	0.005, 0.03	0.005, 0.03	0.005, 0.03	0.005, 0.03	0.005, 0.03
G1	Blinded, GSK	DMSO	50	0.2%	0.05, 0.1	DMSO	50	0.05, 0.1	0.05, 0.1	0.05, 0.1	0.05, 0.1	0.05, 0.1
G2	Blinded, GSK	DMSO	500	0.1%	0.25, 0.5	DMSO	500	0.25, 0.5	0.25, 0.5	0.25, 0.5	0.25, 0.5	0.25, 0.5
G3	Blinded, GSK	DMSO	500	0.02%	0.025, 0.05, 0.1	DMSO	50	0.025, 0.05, 0.1	0.025, 0.05, 0.1	0.025, 0.05, 0.1	0.025, 0.05, 0.1	0.025, 0.05, 0.1

R1	Blinded, Roche	Culture medium	10	-	0.125, 0.5	KHH	10	0.125, 0.5
R2	Blinded, Roche	Culture medium	100	-	0.5, 1	KHH	100	0.5, 1
R3	Blinded, Roche	Milli-Q + 11.4% v/v 1 M NaOH	100	0.1%	0.025, 0.1	Milli-Q + 11.4% v/v 1 M NaOH	100	0.025, 0.1
P1	Blinded, Pfizer	DMSO	500	0.2%	0.25, 1	DMSO	500	500, 1000
P2	Blinded, Pfizer	DMSO	500	0.1%	0.25, 0.5	DMSO	500	0.25, 0.5

CDDP, cisplatin; Conc., concentration; CSA, cyclosporin A; D-F12, Dulbecco's modified Eagle's medium and nutrient mixture F-12; DMSO, dimethyl sulfoxide; HBSS, Hank's balanced salt solution; KHH, Krebs-Henseleit buffer supplemented with HEPES (10 mM); PTEC CM, proximal tubule epithelial cell complete medium; PTEC SFM, proximal tubule epithelial cell serum-free medium; TBR, tobramycin; TNV, tenofovir

Table S2 Overview of available *in vitro* and *in vivo* data of blinded nephrotoxicants

Compound	Previous studies	Reference
G1	<p>Oral dose studies in monkeys (14 days): kidney weights and kidney to body ratios for males dosed at 1500 mg G1 R-enantiomer/kg/day where higher than those for the controls. Especially one animal was particularly high. Mottled and irregular surface was observed in one male kidney; and pallor in one female kidney at 1500 mg/kg/day R-enantiomer. Renal tubular nephritis occurred in two males given 1500 mg/kg/day of G1 R-enantiomer. C_{\max} of R-enantiomer at highest dose 6.9 μM (2.5 $\mu\text{g/mL}$).</p> <p>Nephroscreen flag: 50 μM</p>	Sponsor
G2	<p>Oral toxicity study in rats (1 month): one male and one female died at day 3, two males on day 4. Highest dose of 200 mg/kg/day was lowered to 150 mg/kg/day on day 4. Tubular necrosis associated with tubular regeneration and interstitial mononuclear cell infiltration was noted in three males at the highest dose. One additional animal exhibited serum chemistry changes associated with renal tubular lesions. Thirteen male rats and 17 female rats with similar dosing regimen survived to the end of the study. C_{\max} at highest dose in male was 18 $\mu\text{g/mL}$ (45 μM) and in female 41 $\mu\text{g/mL}$ (103 μM).</p> <p>Oral toxicity study in Beagle dogs (1 month): dose at 0.5, 3 or 20 mg/kg/day. Both kidneys at highest dose in males were reduced in size, correlated with reduced weight and tubular degeneration/regeneration (in 2 out of 3 males) and glomerulosclerosis. C_{\max} at highest dose in male was 14 $\mu\text{g/mL}$ (35 μM) and in female 15 $\mu\text{g/mL}$ (37.5 μM).</p> <p>Nephroscreen flag: MRP interaction, 250 μM</p>	Sponsor

G3

Oral gavage dose range-finding study in rats (7 days): inclusions in epithelial cells lining collecting ducts of papilla, crystals in the collecting ducts. Some animals presented with basophilia of tubules in the cortex and degeneration/regeneration of the medullary epithelium. Some animals (both male and female) presented with tubule dilation. Increases in plasma creatinine concentrations (mild to moderate). C_{\max} in male was 2 µg/mL (4.6 µM) and in female 4.5 µg/mL (10.4 µM).

Sponsor

Oral study in minipig (7 days): 100 or 300 mg/kg/day, once daily; or 1000mg/kg/day up to 3 days (severe clinical signs – animals were killed). Pallor and enlargement of the kidney; dilated tubules, interstitial inflammation/fibrosis in all areas of the kidney, crystals in lumen consistent with obstructive nephropathy. Increased plasma levels of urea and sCr. C_{\max} in male was 19.4 µg/mL (45 µM) and in female 3.5 µg/mL (8 µM).

Nephroscreen flag: 50 µM

R1

[1], [2]

Colistin (antibiotic)

Colistin undergoes extensive renal tubular reabsorption (up to 80%) upon glomerular filtration. Most of the filtered colistin is retained in the body. Clinical use leads to plasma concentrations of approximately 2 mg/L (approx. 1.7 µM) following i.v. administration of 250 mg/day. Colistin increases tubular epithelial barrier permeability and causes cell swelling, cell lysis, resulting in acute tubular necrosis. Conventional 2D human primary PTECs were susceptible for toxicity at 100 µM.

Nephroscreen flag: 250 µM/ 62 µM

R2

[3]

Cefepime

*(cephalosporin
antibiotic)*

Cefepime accumulates in PTECs where it can form crystals and induces necrotic cell death. It is clinically administered at high doses (up to 9 g/day, parenterally). Conventional 2D human primary PTECs were susceptible to high concentrations (500-1000 µM) and after prolonged treatment (several days).

Nephroscreen flag: 500 µM/ 62 µM

R3	Zoledronat (zoledronic acid) is excreted unmetabolized via the kidney. It is administered clinically at doses of 4 mg/day, i.v. . A relationship exists between peak levels of zoledronic acid in the blood and renal toxicity since renal damage declines when the dose is reduced. <i>In vitro</i> toxicity documented in a rat renal cell line (NRK- 52E at 100 µM), human primary cells (100 µM) and in a renal cell line (HK-2 at 1000 µM).			[4]-[7]
Zoledronat (bisphosphonate drug class)				
P1	Vasopressin antagonist	1a	Nephroscreen flag: 25 µM	Sponsor
	<i>In vitro</i> toxicity found in human liver cell line THLE at > 300 µM. In a rat study (2 weeks), females showed tubular degeneration/regeneration at highest dose (125 mg/kg/day). The renal cortex and the outer and inner stripes of the outer medulla were affected. This finding was most suggestive of damage to the proximal tubules but also, to a lesser extent, distal nephron segments, and correlated with slight increases in urinary protein. Proteinuria was most likely the result of decreased resorption of low molecular weight proteins from the proximal convoluted tubule, secondary to tubular damage. C _{max} was (both in male and female) 14.6 µM.			
P2	Prostaglandin receptor antagonist	E2	Nephroscreen flag: 250 µM	Sponsor
	<i>In vitro</i> toxicity found in human liver cell line THLE at > 262 µM. In a rat study (2 weeks), male and female rats were treated with mid dose (50 mg/kg) or high dose (250 mg/kg). In female rats, 1 out of 10 rats at mid dose and all rats at high dose showed renal tubular degeneration in the pars recta, while males were unaffected. Levels of KIM-1 in tissue and urine were elevated, while levels of BUN and sCr remained normal. Full recovery was observed on cessation of dosing. It is noteworthy that AUC in male rats at 50 and 250 mg/kg on day 14 was 6 and 8 times lower, respectively, than those achieved in females at the same dose level based on the free unbound drug fraction. C _{max} was (both in male and female) 1.5 mM, unbound fraction was 55.8 µM.			
	Nephroscreen flag: 250 µM			

2D, two-dimensional, AUC, area under the curve; BUN, blood urea nitrogen; C_{max}, maximum serum concentration; HK-2, human kidney-2; KIM-1, kidney injury molecule-1; NRK- 52E, normal rat kidney-52E; PTEC, proximal tubule epithelial cell; sCr, serum creatinine; THLE, transformed human liver epithelial

References

- [1] Ordooei Javan A, Shokouhi S, Sahraei Z 2015. A review on colistin nephrotoxicity. *Eur J Clin Pharmacol* 71(7):801-810.
- [2] Plachouras D, Karvanen M, Friberg LE, Papadomichelakis E, Antoniadou A, Tsangaris I, Karaiskos I, Poulakou G, Kontopidou F, Armaganidis A, Cars O, Giamarellou H 2009. Population pharmacokinetic analysis of colistin methanesulfonate and colistin after intravenous administration in critically ill patients with infections caused by gram-negative bacteria. *Antimicrob Agents Chemother* 53(8):3430-3436.
- [3] Elsayed MG, Elkomy AA, Gaballah MS, Elbadawy M 2014. Nephrotoxicity of cefepime: A new cephalosporin antibiotic in rats. *J Pharmacol Pharmacother* 5(1):33-38.
- [4] Markowitz GS, Fine PL, Stack JJ, Kunis CL, Radhakrishnan J, Palecki W, Park J, Nasr SH, Hoh S, Siegel DS, D'Agati VD 2003. Toxic acute tubular necrosis following treatment with zoledronate (Zometa). *Kidney Int* 64(1):281-289.
- [5] Luhe A, Kunkele KP, Haiker M, Schad K, Zihlmann C, Bauss F, Suter L, Pfister T 2008. Preclinical evidence for nitrogen-containing bisphosphonate inhibition of farnesyl diphosphate (FPP) synthase in the kidney: implications for renal safety. *Toxicol In vitro* 22(4):899-909.
- [6] Bergner R, Siegrist B, Gretz N, Pohlmeyer-Esch G, Kranzlin B 2015. Nephrotoxicity of ibandronate and zoledronate in Wistar rats with normal renal function and after unilateral nephrectomy. *Pharmacol Res* 99:16-22.
- [7] Verhulst A, Sun S, McKenna CE, D'Haese PC 2015. Endocytotic uptake of zoledronic acid by tubular cells may explain its renal effects in cancer patients receiving high doses of the compound. *PLoS One* 10(3):e0121861.

