



**Universiteit  
Leiden**  
The Netherlands

## **A perfused human blood-brain barrier on-a-chip for high-throughput assessment of barrier function and antibody transport**

Wevers, N.R.; Kasi, D.G.; Gray, T.; Wilschut, K.J.; Smith, B.; Vught, R. van; ... ; Obermeier, B.

### **Citation**

Wevers, N. R., Kasi, D. G., Gray, T., Wilschut, K. J., Smith, B., Vught, R. van, ... Obermeier, B. (2018). A perfused human blood-brain barrier on-a-chip for high-throughput assessment of barrier function and antibody transport. *Fluids And Barriers Of The Cns*, 15(1).  
doi:10.1186/s12987-018-0108-3

Version: Publisher's Version

License: [Creative Commons CC BY 4.0 license](#)

Downloaded from: <https://hdl.handle.net/1887/4105263>

**Note:** To cite this publication please use the final published version (if applicable).

RESEARCH

Open Access



# A perfused human blood–brain barrier on-a-chip for high-throughput assessment of barrier function and antibody transport

Nienke R. Wevers<sup>1,2\*</sup>, Dhanesh G. Kasi<sup>1</sup>, Taylor Gray<sup>3</sup>, Karlijn J. Wilschut<sup>1</sup>, Benjamin Smith<sup>3</sup>, Remko van Vught<sup>1</sup>, Fumitaka Shimizu<sup>4</sup>, Yasuteru Sano<sup>4</sup>, Takashi Kanda<sup>4</sup>, Graham Marsh<sup>3</sup>, Sebastiaan J. Trietsch<sup>1</sup>, Paul Vulto<sup>1</sup>, Henriëtte L. Lanz<sup>1†</sup> and Birgit Obermeier<sup>3†</sup>

## Abstract

**Background:** Receptor-mediated transcytosis is one of the major routes for drug delivery of large molecules into the brain. The aim of this study was to develop a novel model of the human blood–brain barrier (BBB) in a high-throughput microfluidic device. This model can be used to assess passage of large biopharmaceuticals, such as therapeutic antibodies, across the BBB.

**Methods:** The model comprises human cell lines of brain endothelial cells, astrocytes, and pericytes in a two-lane or three-lane microfluidic platform that harbors 96 or 40 chips, respectively, in a 384-well plate format. In each chip, a perfused vessel of brain endothelial cells was grown against an extracellular matrix gel, which was patterned by means of surface tension techniques. Astrocytes and pericytes were added on the other side of the gel to complete the BBB on-a-chip model. Barrier function of the model was studied using fluorescent barrier integrity assays. To test antibody transcytosis, the lumen of the model's endothelial vessel was perfused with an anti-transferrin receptor antibody or with a control antibody. The levels of antibody that penetrated to the basal compartment were quantified using a mesoscale discovery assay.

**Results:** The perfused BBB on-a-chip model shows presence of adherens and tight junctions and severely limits the passage of a 20 kDa FITC-dextran dye. Penetration of the antibody targeting the human transferrin receptor (MEM-189) was markedly higher than penetration of the control antibody (apparent permeability of  $2.9 \times 10^{-5}$  versus  $1.6 \times 10^{-5}$  cm/min, respectively).

**Conclusions:** We demonstrate successful integration of a human BBB microfluidic model in a high-throughput plate-based format that can be used for drug screening purposes. This in vitro model shows sufficient barrier function to study the passage of large molecules and is sensitive to differences in antibody penetration, which could support discovery and engineering of BBB-shuttle technologies.

**Keywords:** Blood–brain barrier, Microfluidics, Organ-on-a-chip, BBB, Antibody transcytosis

\*Correspondence: n.wevers@mimetas.com

†Henriëtte L. Lanz and Birgit Obermeier contributed equally to this work

<sup>1</sup> Mimetas BV, J.H. Oortweg 19, 2333 CH Leiden, The Netherlands

Full list of author information is available at the end of the article



## Background

The blood–brain barrier (BBB) ensures a homeostatic environment for the central nervous system (CNS) and is essential for healthy brain functioning. The BBB comprises specialized endothelial cells and supporting cells, such as astrocytes and pericytes. Due to a combination of specific transport mechanisms and the presence of adherens junctions and tight junctions, the BBB controls passage of compounds into the brain [1–5]. This way, the BBB protects the brain from many harmful substances that circulate in the blood. However, the BBB's barrier properties also complicate the treatment of CNS disorders, as many small- and large-molecule pharmaceuticals are restricted from entering the brain in quantities that are large enough to elicit a therapeutic response [6]. It is therefore necessary to develop improved drug delivery strategies that enable efficient delivery of biopharmaceuticals to the brain.

The BBB employs specialized transporter systems to allow essential nutrients to enter the brain. The transport system that is most attractive to deliver large-molecule drugs into the brain is receptor-mediated transcytosis (RMT). In RMT, a ligand (or antibody) binds a receptor on the luminal surface of a brain endothelial cell, after which it undergoes internalization via endocytosis and is trafficked to the abluminal side, where it can be released and gain access to the brain parenchyma. Harnessing this process for therapeutic drug delivery is compelling, as it could allow for selective transport into the CNS in a non-invasive manner, without disruption of the BBB [7, 8]. Several studies have demonstrated increased CNS exposure to therapeutic antibodies by combining them with RMT targeting antibodies against the transferrin receptor, insulin receptor, low-density lipoprotein receptor-related proteins 1 and 2, and the large neutral amino acid transporter 1 [9–11]. However, challenges exist in optimizing antibody properties (such as affinity, valency, bispecific format, and Fc receptor engagement) to effectively and safely traffic across the brain endothelium [12–16]. Improved *in vitro* models that enable further research of the cellular and molecular mechanisms underlying transcytosis at the BBB are needed to improve these CNS drug delivery technologies [14, 17, 18].

While *in vivo* models can be used to study an intact BBB in its physiological environment, the complexity involved in deciphering whole-organism drug distribution and the lower throughput of these studies limits their use in screening for BBB-penetrant antibodies. For this reason, *in vivo* research in the field is complemented by simpler and faster *in vitro* models, such as the Transwell method [19–22] and several on-a-chip systems [23–29]. Although the field of *in vitro* BBB modelling has tremendously progressed in recent years, there is a need

for a model that combines fast, high-throughput readouts with physiologically relevant conditions, such as flow, co-culture, and the absence of artificial membranes.

In this manuscript, we show the development of an *in vitro* model of the human BBB in a high-throughput microfluidic platform. The platform allows patterning of extracellular matrix gel by means of surface tension. A blood vessel is grown adjacent to that gel and a third channel is used to insert astrocytes and pericytes. The system is free of artificial membranes, accommodates fluid flow through the blood vessels, and allows fluid-phase sampling of molecules that penetrate the endothelial and matrix layers. Using two different antibodies, we show that the model is sensitive to differences in antibody penetration of brain endothelial cells. The model may support further discovery of antibody BBB-shuttle technologies.

## Methods

### Cell culture

Cell lines of brain endothelial cells, pericytes, and astrocytes were provided by Yamaguchi University, Japan, and originate from the following human primary cell sources: human brain microvascular endothelial cells (TY10 cell line) were isolated from normal brain tissue from a patient with meningioma. Human brain pericytes (hBPCT cell line) were derived from brain tissue of a patient that died from a heart attack. Human astrocytes (hAst cell line) were generated from human primary astrocytes distributed by Lonza (Basel, Switzerland). All three cell types were immortalized with retroviral vectors harboring a SV40 large T antigen gene that is engineered to drive proliferation at 33 °C [30–34] and have been used to model the BBB in previous studies [35–38]. Cells were cultured at 33 °C, 5% CO<sub>2</sub> to allow optimal cell expansion in T75 flasks (734-2705, Corning, NY, USA), which were pre-coated with 50 µg/mL collagen-I (Cultrex 3D collagen-I Rat Tail, 5 mg/mL, 3447-020-01, AMS-bio, Abingdon, UK) in 1% acetic acid (A6283, Sigma, St. Louis, MO, USA) in water. TY10 cells were used between passage 17–25 and cultured in ScienCell endothelial cell medium (#1001, ScienCell, Carlsbad, CA, USA). The hAst cells were used between passage 7–12 and cultured in ScienCell astrocyte medium (#1018, ScienCell). The hBPCT cells were used between passage 14–25 and cultured in ScienCell pericyte medium (#1012, ScienCell). Cells were routinely tested for mycoplasma contamination and found negative.

### Culture of TY10 microvessels in the two-lane OrganoPlate

Two-lane OrganoPlates (Mimetas BV, the Netherlands) with 400 µm × 220 µm (w × h) channels were employed. Phaseguides had dimensions of 100 µm × 55 µm (w × h)

and 2  $\mu\text{L}$  of gel composed of 4 mg/mL collagen-I (Cultrex 3D collagen-I Rat Tail, 5 mg/mL, 3447-020-01, AMSbio), 100 mM HEPES (15630-122, Thermo Fisher, Waltham, MA, USA) and 3.7 mg/mL  $\text{NaHCO}_3$  (S5761, Sigma) was dispensed in the gel inlet and the OrganoPlate was incubated for 15 min at 33 °C. After plate incubation, 25  $\mu\text{L}$  of PBS was added to the gel inlet to prevent the gel from drying out. A TY10 cell suspension of  $1.5 \times 10^7$  cells/mL was prepared and 2  $\mu\text{L}$  was seeded in the medium inlet. 50  $\mu\text{L}$  of medium was added to the medium inlet and PBS was aspirated from the gel inlet. The plate was incubated on the side for 3 h in the incubator to allow the cells to sediment against the collagen-I gel and attach. After incubation, 50  $\mu\text{L}$  of medium was added to the medium outlet. The OrganoPlate was placed on an interval rocker switching between a +7° and −7° inclination every 8 min (Mimetas Rocker Mini, Mimetas BV), allowing bidirectional flow. Cells were cultured at 33 °C (and 5%  $\text{CO}_2$ ) to allow full cell coverage of the ECM gel. Medium was refreshed every 2–3 days. A schematic representation of all steps is shown in Additional file 1. The following media were used to assess their influence on barrier function: ScienCell endothelial cell medium (#1001, ScienCell), Cell Biologics endothelial cell medium (#H1168, Cell Biologics, Chicago, IL, USA), MV2 medium (C-22121, Bioconnect, Huissen, the Netherlands), and EBM-2 medium (cc-3156, Lonza).

#### BBB co-culture in the three-lane OrganoPlate

OrganoPlate BBB co-culture was performed using three-lane OrganoPlates with  $400 \mu\text{m} \times 220 \mu\text{m}$  ( $w \times h$ ) channels (Mimetas BV). Phaseguides had dimensions of  $100 \mu\text{m} \times 55 \mu\text{m}$  ( $w \times h$ ). To establish a BBB co-culture, a collagen-I gel was dispensed in the gel inlet of the chips and filled the middle channel. TY10 cells were seeded and cultured in the top channel as described in the previous section. After 3 days, 2  $\mu\text{L}$  of a  $7 \times 10^6$  cells/mL cell suspension of hAst and hBPCT cells (1:3 ratio) was seeded in the bottom channel. The plate was incubated on the side for 1.5 h to allow the hAst and hBPCT cells to attach to the collagen-I gel. After incubation, fresh ScienCell astrocyte medium is added to the inlet and outlets of the top channel (50  $\mu\text{L}$  in each), after which perfusion is reinstated by placing the plate on the rocker platform. During the entire culture period, only the top channel, which contains the TY10 microvessel, is perfused to allow optimal endothelial barrier strength. Medium was refreshed every 2–3 days. Assays were performed at day 7. A schematic representation of all steps is shown in Additional file 2.

For the images shown in Fig. 3f, g, hBPCT cells were labeled with Calcein red<sup>TM</sup> AM (21900, AAT Bioquest, Sunnyvale, CA, USA) and hAst cells were labeled with

green-fluorescent calcein-AM (C3099, Thermo Fisher) before seeding in the OrganoPlate.

#### Immunocytochemistry

Cultures in the OrganoPlate were fixed with 100% methanol (−20 °C, 494437, Sigma) and permeabilized with 0.3% Triton X-100 (T8787, Sigma). Cells were incubated with blocking solution (2% FCS, 2% bovine serum albumin (BSA, A2153, Sigma), 0.1% Tween-20 (P9416, Sigma)) for 45 min. Primary antibody was incubated for 1–2 h, after which secondary antibody was incubated for 30 min. The following antibodies were used: anti-claudin-5 (35-2500, Thermo Fisher), anti-VE-cadherin (ab33168, Abcam), anti-PECAM-1 (M0823, Dako), goat anti-rabbit AlexaFluor 488 (A11008, Thermo Fisher), goat anti-rabbit AlexaFluor 555 (A21428, Thermo Fisher), goat anti-mouse AlexaFluor 488 (A11001, Thermo Fisher), goat anti-mouse AlexaFluor 555 (A21422, Thermo Fisher), and donkey anti-mouse AlexaFluor 647 (A31571, Thermo Fisher). Nuclei were stained using Hoechst (H3570, Thermo Fisher). All steps were performed at room temperature (RT). Cells were imaged with ImageXpress Micro XLS and Micro XLS-C HCI Systems (Molecular Devices, San Jose, CA, USA).

TY10 cells grown in a collagen-I coated 24-well glass bottom plate (P24-0-N, Cellvis, Mountain View, CA, USA) were fixed with 4% paraformaldehyde (50-980-487, Thermo Fisher) for 10 min and incubated with primary antibody against the human transferrin receptor (A11130, Thermo Fisher) for 3 h, followed by 1 h incubation with goat anti-mouse Alexa Fluor 488 (A11001, Thermo Fisher). Nuclei were stained with Hoechst and cells were imaged using a Zeiss LSM 710 Confocal Microscope (Zeiss, Oberkochen, Germany).

#### Barrier integrity assay

Chips were washed with culture medium (25  $\mu\text{L}$  on all inlets and outlets,  $1 \times 5$  min) to ensure proper flow profiles during the subsequent barrier integrity assay. Next, all medium was aspirated from the chips and 20  $\mu\text{L}$  of medium without fluorescent compound was added to the basal side of the chips (for the two-lane OrganoPlate this is the gel inlet, for the three-lane OrganoPlate these are the gel inlets and outlets and bottom medium inlets and outlets). Medium containing 0.1 mg/mL FITC-dextran (20 kDa, FD20S, Sigma) was added to the top channel, which contained the TY10 microvessel (40  $\mu\text{L}$  on inlet, 30  $\mu\text{L}$  on outlet) and image acquisition was started. Leakage of the fluorescent molecule from the lumen of the microvessel into the adjacent gel channel was automatically imaged using an ImageXpress XLS Micro HCI system (molecular devices). The ratio between the fluorescent signal in the basal and apical region of the tube

was analyzed using Fiji [39]. Graphs were plotted using GraphPad Prism 6 (GraphPad Software, San Diego, CA, USA).

#### Analysis of cell surface binding of anti-hTfR MEM-189

TY10 cells were cultured to confluency and lifted with accutase for 1 h.  $2.5 \times 10^5$  cells/mL were mixed with antibody in PBS with 0.5% BSA and incubated for 1 h at 4 °C. Cells were washed 3× and incubated with 3 µg/mL PE-goat anti-mIgG (115-116-146, Jackson ImmunoResearch, Cambridgeshire, UK) for 1 h at 4 °C, then washed and fixed with 1% paraformaldehyde for 10 min at RT. Cell fluorescence intensity was analyzed by flow cytometry (FACSCalibur, BD Biosciences, Franklin Lakes, NJ, USA) and data was analyzed in FlowJo v10 (FlowJo LLC, Ashland, OR, USA) and GraphPad Prism software.

#### Antibody transcytosis assay

Anti-human transferrin receptor mouse monoclonal antibody MEM-189 mIgG1 (MA1-21562, Thermo Fisher, 10 × 0.1 mg) was dialyzed into pyrogen free PBS to remove the azide in the supplied storage solution. For negative control, an anti-hen egg lysozyme (anti-HEL) antibody (F10.6.6, Genbank AF110316 VH and AY277254.1 VL) was expressed as a mouse IgG1 antibody in CHO cells and purified by recombinant protein A Sepharose (GE) affinity chromatography and size exclusion chromatography (superdex 200).

Chips were washed once with medium to ensure proper flow profiles (as described for the barrier integrity assay). Next, medium was aspirated from the chips, after which 20 µL of fresh medium without antibody was added to the gel inlets and outlets and the bottom medium inlets and outlets. A total of 70 µL of antibody dilution (1.25 µM in medium) was added to the top channel and the OrganoPlate was incubated on the rocker platform in the incubator for 1 h after which basal samples were taken, which consisted of the full contents of the gel inlet and outlet and bottom medium inlet and outlet.

#### Analysis of antibody contents in basal samples

A Meso Scale Discovery (MSD) platform based quantitative immunoassay was used to determine the concentration of antibodies in basal samples collected from the OrganoPlate. Multi-Array 96-well MSD high-binding plates (L15XB-3/L11XB-3, Meso Scale Discovery, Rockville, MD, USA) were coated overnight at 4 °C with 5 µg/mL capturing agent AffiniPure goat anti-mouse IgG, Fcγ fragment specific (115-005-071, Jackson ImmunoResearch). Plates were blocked (1% BSA in PBS) for 2 h, after which 100 µL/well standards (generated from 1.25 µM antibody stock dilutions) or samples were added to the plate and incubated for 1.5–2 h. Plates were

incubated with 100 µL/well 0.25 µg/mL primary detection agent Biotin-SP-conjugated AffiniPure goat-anti-mouse IgG (115-005-071, Jackson ImmunoResearch) for 45 min, followed by 100 µL/well 0.75 µg/mL secondary detection agent MSD Sulfo-TAG StreptAvidin (R32AD-1, Meso Scale Discovery) for 30 min. Between each step, the plates were washed 4× with wash buffer (PBST: 1× PBS with 0.05% Tween-20, 28352, Thermo Fisher). Immediately before MSD read, 2× Read Buffer (MSD Read Buffer T (4×), R92TD-2, diluted 1:2 in diH<sub>2</sub>O) was added to the plates. Plates were read on a MSD QuickPlex SQ 120 in automatic read mode, then a Sigmoidal, 4PL,  $X = \log(\text{concentration})$  interpolation was used to determine antibody concentration in GraphPad Prism 7.02 (GraphPad Software). The apparent permeability ( $P_{\text{app}}$ ) of both antibodies was determined using the following formula:  $P_{\text{app}} = [\Delta C_{\text{receiver}} / \Delta t] \times [V_{\text{receiver}} / (A_{\text{barrier}} \times C_{\text{donor, initial}})]$ , in which  $\Delta C_{\text{receiver}} / \Delta t$  is the change in antibody concentration in the receiving compartment over time,  $V_{\text{receiver}}$  is the volume of the receiving compartment,  $A_{\text{barrier}}$  is the surface area of the barrier, and  $C_{\text{donor, initial}}$  is the initial antibody concentration of the donor compartment.

## Results


#### Perfused microvessels of TY10 brain endothelial cells in the two-lane OrganoPlate

The two-lane OrganoPlate (see Fig. 1a, b) allows parallel culture of 96 miniaturized tissues on microfluidic tissue chips [40–42]. In each chip, a microvessel of brain endothelial cells was grown under perfusion against an extracellular matrix gel (see Fig. 1c–e and Additional file 1). A 3D reconstruction of a TY10 microvessel stained with ActinRed is shown in Fig. 1e and demonstrates that a vessel-like structure of endothelial cells has formed. Immunostaining was performed to assess the expression of key BBB adherens junction and tight junction proteins (see Fig. 1f–h). TY10 microvessels showed expression of vascular endothelial cadherin (VE-cadherin), platelet endothelial cell adhesion molecule 1 (PECAM-1), and claudin-5 at the cell–cell contacts as expected. The expression and interendothelial localization of these markers are indicative of barrier formation [2, 43].

#### Assessment of barrier function in TY10 brain endothelial microvessels

The barrier function of TY10 microvessels grown in the OrganoPlate was assessed using a fluorescent barrier integrity assay. TY10 microvessels were perfused with a 20 kDa FITC-dextran dye, which has a hydrodynamic radius slightly smaller than a folded IgG antibody (3 nm [44] versus 5–6 nm [45], respectively). Leakage of the dye from the microvessel into the adjacent gel channel



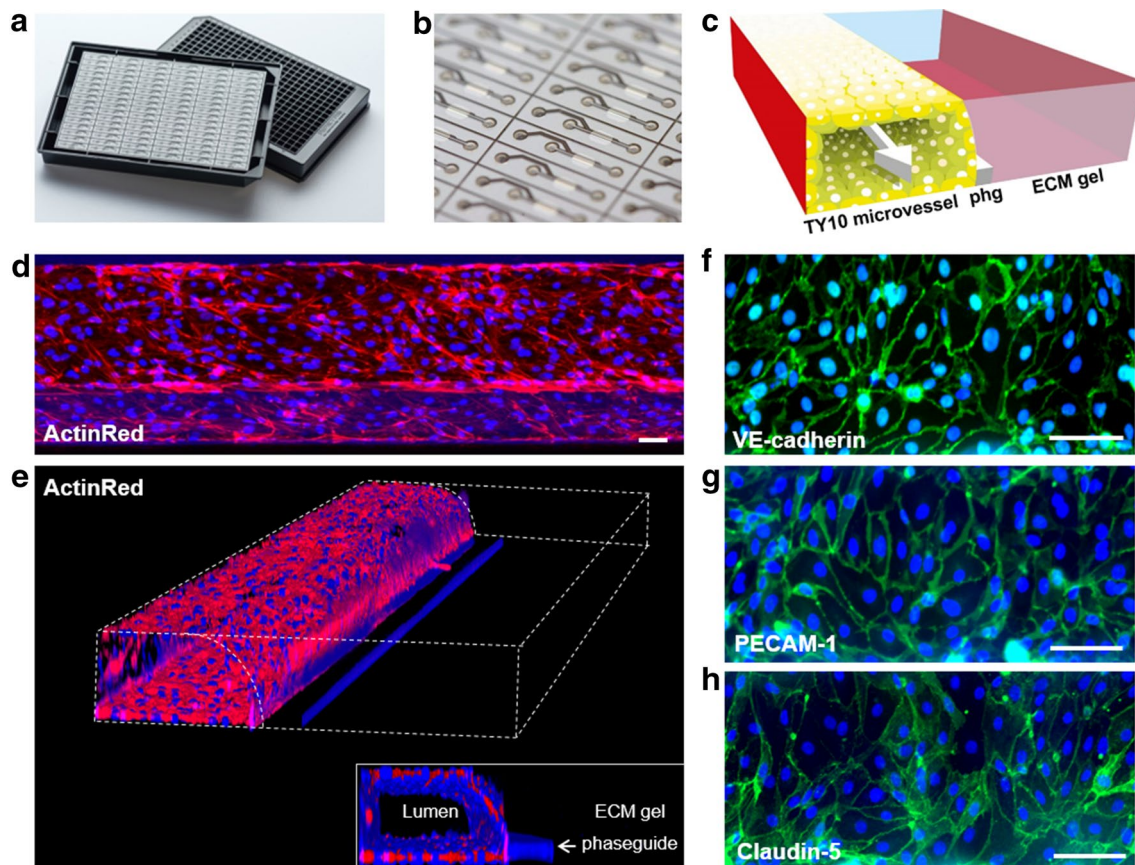


**Fig. 1** Perfused microvessels of TY10 brain endothelial cells in the two-lane OrganoPlate. **a** Picture shows the two-lane OrganoPlate. The plate combines a 384-well microtiter plate on the top with microfluidic channels on the bottom that make up 96 tissue culture chips. **b** Zoom-in of the bottom of the two-lane OrganoPlate, showing the tissue culture chips that consist of two channels: a gel channel and a medium channel. **c** 3D artist impression of the center of a chip. ECM gel is added to the gel channel and a phaseguide (phg) prevents it from flowing into the adjacent medium channel. After gelation of the ECM gel, TY10 cells are added to the medium channel and a TY10 microvessel forms. The microvessel has a lumen at its apical side that is perfused. **d** Maximum projection image of a TY10 microvessel stained with ActinRed. Scale bar is 50  $\mu$ m. **e** 3D reconstruction of a confocal z-stack showing a perfused TY10 microvessel. Inlay shows a vertical cross section, depicting the lumen, the phaseguide, and the ECM gel channel. **f-h** Immunostaining of a TY10 microvessel for VE-cadherin, PECAM-1, and claudin-5. Scale bar is 50  $\mu$ m

was assessed by acquisition of fluorescent images over time (see Fig. 2a). In a chip containing a leak-tight TY10 microvessel, all fluorescent dye is retained

within the vessel (see Fig. 2b, top image). In a cell-free control chip, the fluorescent dye can freely diffuse into the adjacent gel channel (see Fig. 2b, middle image).

(See figure on next page.)

**Fig. 2** Assessment of barrier function in TY10 brain endothelial microvessels. **a** Schematic representation of the barrier integrity assay. A perfused TY10 microvessel is grown against an ECM gel. A fluorescent dye is added to the medium inlets and outlets and perfused through the lumen of the microvessel. In case of a leaktight vessel, all dye is retained in the vessel. In case of a leaky vessel, the dye leaks into the adjacent gel channel. **b** Fluorescent dye distribution for a TY10 microvessel (top image) and a cell-free control (middle image) during a barrier integrity assay (FITC-dextran, 20 kDa,  $t = 30$  min). To quantify leakage of the fluorescent dye in a chip, the fluorescence intensity measured in the medium channel and the gel channel ( $Fluo_{Med}$  and  $Fluo_{Gel}$ , respectively) are measured over time (bottom image). **c** For each condition the ratio between the fluorescence signals measured in the two channels is plotted over time. In case of a leaktight TY10 microvessel, the ratio  $Fluo_{Gel}/Fluo_{Med}$  is constant, resulting in a flat horizontal line. In case of a leaky microvessel or a cell-free control, the fluorescence signal measured in the gel channel increases over time, resulting in an increase in the ratio  $Fluo_{Gel}/Fluo_{Med}$ . **d** TY10 microvessels were grown under perfusion or static conditions for 7 days, after which barrier integrity assays (for FITC-dextran, 20 kDa) were performed.  $n = 6$  for TY10 microvessels,  $n = 2$  for cell-free controls. Error bars show standard deviation of the mean. **e** Assessment of barrier function (for FITC-dextran, 20 kDa) of TY10 microvessels cultured under perfusion for 7 days in various commercially available cell culture media.  $n = 7$  for all conditions. Error bars show standard deviation of the mean

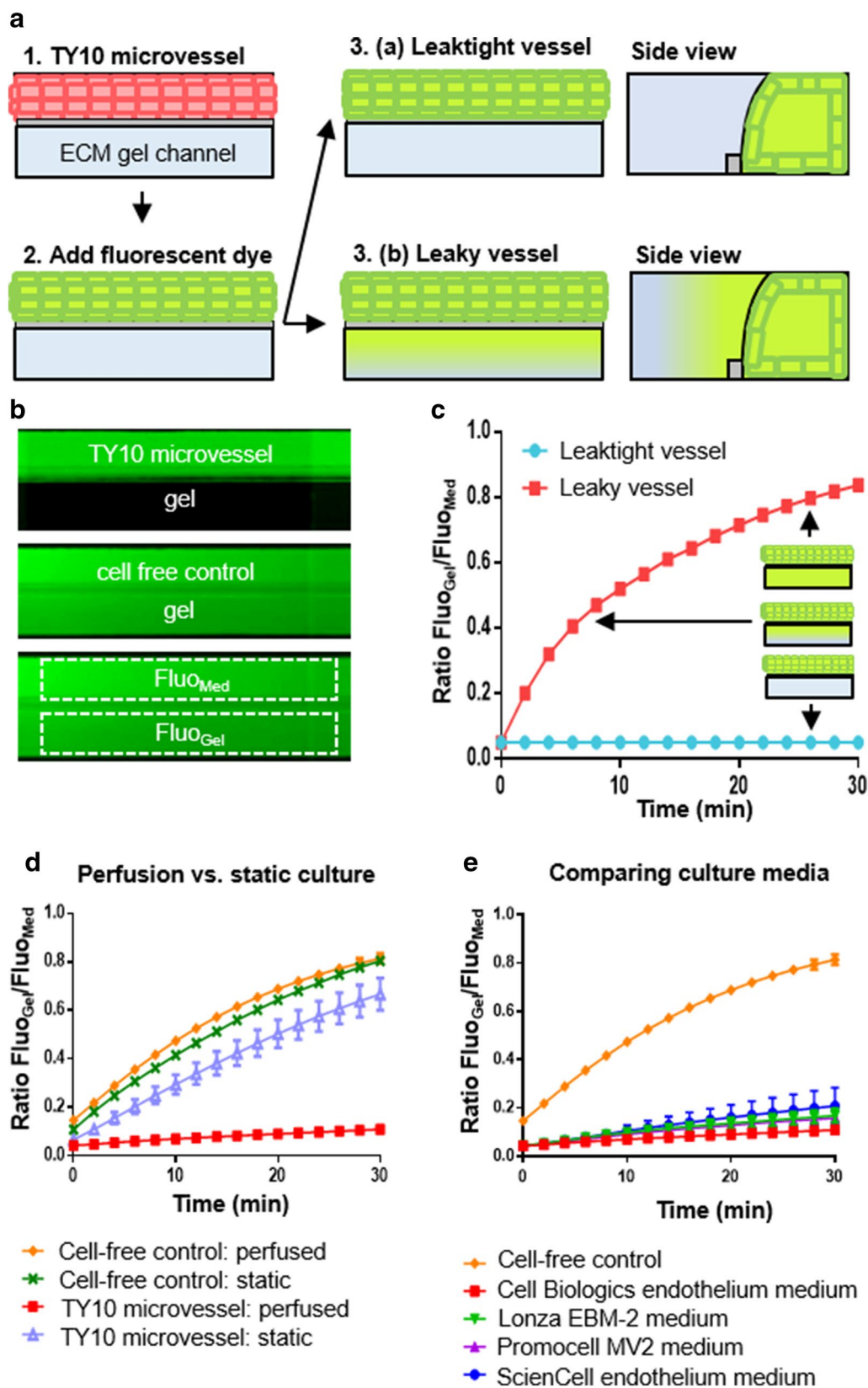

To quantify leakage of the fluorescent dye, areas were selected automatically for each chip to compare over time (see Fig. 2b, bottom image). The ratio of the fluorescence signal measured in the gel channel and the medium channel is plotted over time as an average of all chips within a condition (see Fig. 2c). In chips that contain a tight cellular barrier, this ratio remains constant, as all dye is retained in the microvessel that is grown in the medium channel. In leaky vessels and cell-free controls, the ratio increases over time, as dye leaks into the gel channel. We compared barrier integrity of TY10 microvessels cultured under static conditions to microvessels cultured under perfusion (by passive leveling using an interval rocker). TY10 cells cultured under perfusion show much tighter barrier formation (Fig. 2d) resulting from increased proliferation and elevated expression and improved localization of junctional proteins (see Additional file 3). In addition, we assessed four different commercially available endothelial cell culture media and found that for this specific type of brain microvascular endothelial cell, the endothelial cell medium from Cell Biologics was optimal for barrier formation (see Fig. 2e).

#### BBB co-cultures of brain endothelium, astrocytes, and pericytes in the three-lane OrganoPlate

While endothelial cells make up the brain's vasculature, other cell types, such as astrocytes and pericytes, are also part of the BBB and help maintain barrier function. A three-lane OrganoPlate [46] (see Fig. 3a, b) was employed to establish a BBB co-culture of TY10 brain endothelial cells with hAst and hBPCT cells (see Fig. 3c, f, g and Additional file 2). A barrier integrity assay (see Fig. 3h) was performed to determine barrier function in TY10 monocultures (Fig. 3d) and BBB co-cultures (Fig. 3e) at different time points during culture. Both TY10 monocultures and BBB co-cultures successfully retained, to a great extent, a 20 kDa FITC-dextran dye in the endothelial microvessel at days 5, 7, and 9 as is apparent by the marginal increase in the ratio of fluorescence signal measured in the gel channel and the medium channel over time (see Fig. 3i). This barrier integrity assay is useful in identifying optimal culture conditions and timing for which to perform antibody transcytosis assays and showed minimal passive permeability of molecules that are similar in size to antibodies.

#### Assessment of antibody transcytosis across the BBB


Two different antibodies were used to assess their passage across our human in vitro BBB-on-a-chip. The first was MEM-189, an antibody that binds the human

transferrin receptor (hTfR), which is expressed by TY10 endothelial cells (see Additional file 4a) and has been reported to undergo RMT [47]. Flow cytometry analysis showed that anti-hTfR MEM-189 bound to TY10 cells and that binding was not blocked by 25  $\mu\text{g/mL}$  transferrin [ $\text{EC}_{50} = 0.44 \pm 0.09 \text{ nM}$  (–Tf);  $0.5 \pm 0.1 \text{ nM}$  (+Tf)], see Additional file 4b). The second antibody, anti-hen egg lysozyme (anti-HEL), does not bind a target on human cells and served as a negative control [48]. Both antibodies showed the expected molecular weight, heavy-light chain composition, and size in solution for properly folded antibodies. The antibodies were perfused through the lumen of BBB co-cultures in the three-lane OrganoPlate for 1 h, after which samples were collected from the basal side of the chips (see Fig. 4a). Antibody concentrations in basal samples were determined using MesoScale Discovery® MULTI-ARRAY® technology and the apparent permeability ( $P_{\text{app}}$ ) was calculated. Similar basal concentrations of both antibodies were measured in samples taken from chips without a TY10 microvessel, indicating that when no barrier is present, both antibodies diffuse through the gel to an equal extent (see Fig. 4b). However, in BBB co-cultures, passage of antibody MEM-189 was markedly higher than passage of the control antibody ( $2.9 \times 10^{-5}$  versus  $1.6 \times 10^{-5} \text{ cm/min}$ , respectively) (see Fig. 4b).

#### Discussion

This study describes a significant advancement in the development of a novel BBB model that incorporates human brain endothelial cells, astrocytes, and pericytes in a high-throughput microfluidic platform that can be used for screening purposes. The cells used in this model are immortalized and thus their phenotype may differ from cells of the BBB in living organisms. Therefore, results obtained in this model may not directly apply to patients. However, the TY10 brain endothelial cell line used in this study has been shown to express relevant junctional markers and transporters, independent of the passage number [31]. Although primary brain endothelial cells offer an interesting alternative, they are difficult to obtain reliably from human donors and have been described to rapidly dedifferentiate and lose their BBB characteristics after removal from the in vivo environment, resulting in decreased expression of essential BBB modulators and impaired barrier function [49, 50]. In recent years, tremendous progress has been made in the generation of induced pluripotent stem cell (iPSC)-derived brain endothelial cells. These cells have been shown to express relevant junctional proteins and transporters, respond well to cues from supporting cell types



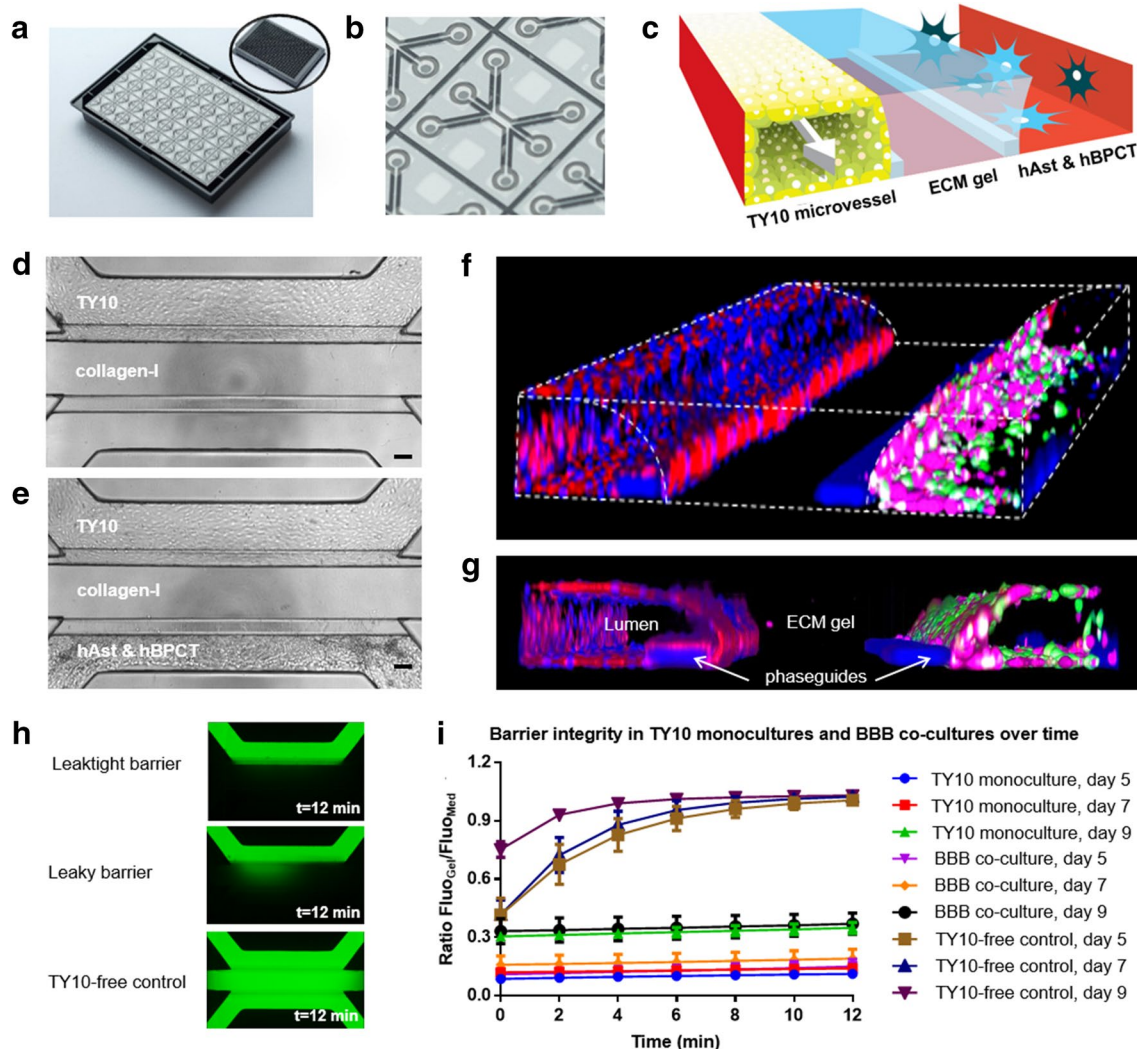


**Fig. 3** BBB co-cultures of brain endothelium, astrocytes, and pericytes in the three-lane OrganoPlate. **a** Picture shows the three-lane OrganoPlate. The plate combines a 384-well microtiter plate on the top with microfluidic channels on the bottom that make up 40 tissue culture chips. **b** Zoom-in of the bottom of the three-lane OrganoPlate, showing one tissue culture chip that consists of three channels. **c** 3D artist impression of the center of a chip. ECM gel is added to the gel channel and two phaseguides (white rims) prevent it from flowing into the adjacent medium channels. After gelation of the ECM gel, TY10 cells are added to one of the medium channels and allowed to form a vessel structure. The microvessel has a lumen at its apical side that is perfused. Next, astrocytes (hAst cells) and pericytes (hBPCT cells) are added to the other medium channel to establish a BBB co-culture. **d, e** Phase contrast images of a TY10 monoculture (top image) and a BBB co-culture (bottom image) in the three-lane OrganoPlate, day 7. Scale bar is 100  $\mu$ m. **f, g** 3D reconstruction of a confocal z-stack showing the organization of the three cell types in a BBB co-culture. The hAst and hBPCT cells were labeled with green-fluorescent calcein-AM and calcein red<sup>TM</sup> AM, respectively, here depicted in the colors green and magenta. Tight junctions of the TY10 microvessel are shown by a PECAM-1 staining (red). **h** Images acquired of chips with a leaktight barrier (top), a leaky barrier (middle), or TY10-free control (gel with hAst and hBPCT, bottom) during a barrier integrity assay (for FITC-dextran, 20 kDa,  $t = 12$  min). **i** Assessment of barrier function (for FITC-dextran, 20 kDa) at different time points for TY10 monocultures and BBB co-cultures cultured under perfusion in the three-lane OrganoPlate.  $n = 3$  for TY10-free controls and TY10 monocultures and  $n = 6$  for BBB co-cultures. Error bars show standard deviation of the mean

such as astrocytes and pericytes, and display in vivo-like barrier properties [51–55]. In addition, a recent report showed that iPSC-derived brain endothelial cells could be used to study antibody transcytosis across the BBB [55]. The combination of our high-throughput,

membrane-free, perfused platform with the ongoing advancements in iPSC-derived cell types may potentially bring about highly relevant BBB models in the future.

In contrast to the standard 2D Transwell approach [19, 56–58], the OrganoPlate supports perfused culture

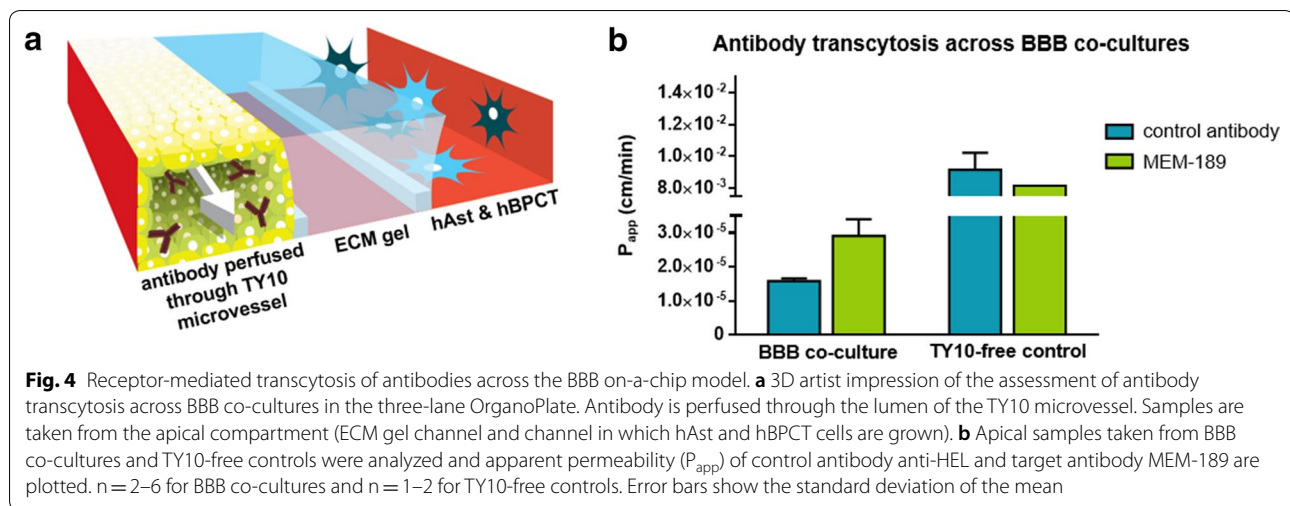

of brain endothelial cells, which was shown to be essential for proper lumen formation, proper expression and localization of junctional proteins, and improved barrier function (see Fig. 2d and Additional file 3). Perfusion is generated by placing the entire plate on a rocker platform, allowing medium to flow from medium inlet to outlet and back, creating a bidirectional flow. Fluid flow is controlled by regulating the inclination angle and the interval with which the rocker platform switches sides. Although direct *in vivo* measurements of shear stress in BBB vessels and venules of different diameters are lacking, it is likely that the shear stress established in our model ( $\sim 1.2$  dyne/cm<sup>2</sup>) is low compared to the shear stress experienced by vessels of similar diameter and curvature *in vivo* [59, 60]. Several reports have described improved barrier function of brain endothelial cells as a result of increased shear stress due to an upregulation in junctional proteins [61–63]. Interestingly, other studies did not find increased expression of junctional proteins, but report that brain endothelial cells exhibit a unique response to shear stress compared to endothelium in different organs. Unlike other endothelial cell types, brain endothelial cells did not elongate or align in response to shear stress, a phenotype that may be associated with the BBB's unique properties [64, 65]. The setup used in this study bypasses the need for pumps and intricate tubing systems, which are often associated with microfluidic culture systems, and improves the user-friendliness and throughput of the method.

Comparison of barrier integrity in different culture platforms is challenging, as the properties of the platform itself often influence the measured outcome and are difficult to properly correct for. Among these properties are the presence or lack of a membrane, the pore size of the membrane, the volumes used for the assay, the

presence of flow, and the nature of the read-out. Standard measures such as a barrier's transelectrical endothelial resistance (TEER) or a compound's  $P_{app}$  can be used to compare results obtained within the same platform, i.e. a Transwell, but cannot directly be compared to results obtained in other culture platforms. The observation that the BBB on-a-chip successfully limits the passage of molecules of similar size as antibodies shows that sufficient barrier function is established to investigate antibody passage. In addition, Trietsch et al. [46] have reported higher sensitivity in detection of compound-induced CaCo-2 barrier disruption in the OrganoPlate compared to Transwell. The increased sensitivity was expected to result from improved maturity of the culture as well as a decreased dead volume and higher surface-to-volume ratio in the OrganoPlate.

Many microfluidic systems employ polydimethylsiloxane (PDMS) because of its optical properties and ease of use in the fabrication process. However, PDMS has several limitations in organ-on-a-chip applications. The intrinsic hydrophobicity of the material impedes cell adhesion and can cause non-specific absorption of proteins and hydrophobic analytes. Although several methods are available to reduce PDMS hydrophobicity and fouling problems, the most successful method requires a complex process that is difficult to incorporate for large scale production [66, 67]. The OrganoPlate employs optical quality (170  $\mu$ m) glass and polymers that are bio-compatible and low compound-absorbing and does not include an artificial membrane. Furthermore, the collagen gel does not significantly restrict passage of antibody into the basal compartment, which was demonstrated by the observation that antibody  $P_{app}$  increased >500-fold when the endothelial cells were omitted from the model (Fig. 4b). Together, the low level passive permeability

of the brain endothelial cell layer, minimal absorption, and the ease of antibody sampling from the basal compartment could support sensitive and high-throughput screening of antibody transcytosis in this BBB on-a-chip.

The human transferrin receptor is of special interest for drug targeting to the brain due to its expression on brain endothelial cells and potential to support receptor mediated transcytosis [9, 12, 68, 69]. In this study, we observed an approximately two-fold higher passage of MEM-189, an antibody targeting the hTfR receptor, across our BBB model compared to a control antibody. The  $P_{app}$  value we measured for the non-binding control antibody agrees well with those reported recently in other in vitro human BBB models [55, 70]. The passage of this murine IgG control antibody could in theory result from Fc receptor mediated transport through the neonatal Fc receptor (FcRn), which is expressed at the BBB. However, since murine IgG1 antibodies show very little binding to the human FcRn [71] and FcRn likely does not result in BBB transcytosis [72], the passage of this antibody is most likely the result of paracellular flux or non-receptor mediated endocytic flux (i.e. micropinocytosis). The enhanced permeability of the MEM-189 antibody is consistent with previous reports [47] and could result from active transport mediated by the transferrin receptor. This bivalent antibody has high affinity for hTfR and thus is not optimized for high transport, as has been shown for other anti-TfR BBB shuttles [9, 12, 15]. As a large portion of endocytosed TfR has been shown to remain in the endothelial cells instead of undergoing transcytosis, it is likely that a relatively large quantity of endocytosed MEM-189 remains within the endothelial cells [73]. However, Sade et al. [47] have shown that binding of MEM-189 to the TfR is pH-dependent, which may support transcytosis via release from the TfR in the endosomal environment and could contribute to its observed transcytosis. The assay developed in this study is sensitive to differences in antibody passage across the BBB model and holds potential to be applied in larger screens for discovery of CNS-penetrant antibodies. The addition of neurons to the model may add an extra layer of complexity that brings about possibilities for modelling the full neurovascular unit and studying the passage of molecules across the BBB as well as the effects of these molecules on neuronal function.

## Conclusion

In summary, the model described here is the first perfused BBB-on-a-chip culture system that is compatible with standard laboratory equipment and allows throughput to an extent that is necessary for drug screening. Moreover, co-culture complexity is easily expanded as shown by the addition of pericytes and astrocytes. This

model could support further discovery and engineering of antibody BBB-shuttle technologies.

## Additional files

**Additional file 1.** Endothelial microvessel seeding in the two-lane OrganoPlate. (a) Schematic representation of one chip of a two-lane OrganoPlate. (b) An ECM gel is seeded in the gel channel, after which endothelial cells are seeded in the medium channel. (c) Endothelial cells attach to the ECM gel and perfusion is started by placing the OrganoPlate on a rocker platform. (d) A microvessel of endothelial cells is formed. (e–g) Cross sectional view of steps described in b–d.

**Additional file 2.** BBB co-culture seeding in the three-lane OrganoPlate®. (a) Schematic representation of one chip of a three-lane OrganoPlate. (b) ECM gel is seeded in the middle gel of the chip, after which endothelial cells (TY10) are seeded in the top channel. (c) Endothelial cells attach to the ECM and perfusion is started by placing the plate on a rocking platform. (d) A microvessel of endothelial cells forms in the top channel, against the ECM gel. (e) Astrocytes (hAst) and pericytes (hBPCTs) are seeded in the bottom channel. (f) hAst and hBPCT cells attach and a BBB co-culture is established. (g–k) Cross sectional view of steps described in b–f.

**Additional file 3.** Comparing perfused and static culture of TY10 microvessels. (a, b) Phase contrast images of TY10 microvessels grown in the two-lane OrganoPlate under perfused or static conditions (day 7). Scale bar is 100  $\mu$ m. (c) Microvessels grown under perfused or static conditions were fixed and nuclei were stained with Hoechst. The average number of nuclei was counted in both conditions and normalized to the perfused condition.  $n = 6$ , Student's t-test  $p < 0.05$ . (d–f) Immunofluorescent staining of TY10 microvessels grown under perfusion for adherens and tight junction markers VE-cadherin, claudin-5, and PECAM-1. (g–i) Immunofluorescent staining of TY10 microvessels grown static for adherens and tight junction markers VE-cadherin, claudin-5, and PECAM-1. Scale bar is 100  $\mu$ m.

**Additional file 4.** Characterization of the human transferrin receptor in TY10 endothelial cells. (a) Immunofluorescent staining of the hTfR in TY10 endothelial cells. Scale bar is 50  $\mu$ m. (b) Flow cytometry analysis of cell surface binding of anti-TfR MEM-189 to TY10 endothelial cells in the presence and absence of transferrin (25  $\mu$ g/mL),  $EC_{50} = 0.44 \pm 0.09$  nM (–Tf);  $0.5 \pm 0.1$  nM (+Tf).

## Abbreviations

BBB: blood–brain barrier; CNS: central nervous system; FcRn: neonatal Fc receptor; FCS: fetal calf serum; hAst: human astrocyte cell line; hBPCT: human brain pericyte cell line; HEL: hen egg lysozyme; HEPES: 4-(2-hydroxyethyl)-1-piperazineethanesulfonic acid; hTfR: human transferrin receptor; iPSC: induced pluripotent stem cell; MEM-189: anti-transferrin receptor antibody; PBS: phosphate buffered saline; PDMS: polydimethylsiloxane; PECAM-1: platelet endothelial cell adhesion molecule 1; RMT: receptor-mediated transcytosis; TEER: transendothelial electrical resistance; TY10 cell line: human brain endothelial cell line; VE-cadherin: vascular endothelial cadherin.

## Authors' contributions

Authors NRW, BS, RV, GM, and BO designed the study. Authors NRW, DGK, TG, and BS performed experiments and data analysis. Authors KJW, SJT, PV, HLL, and BO supervised the research. Authors FS, YS, and TK developed the cells and advised on cell culture procedures. Authors NRW, DGK, TG, BS, HLL, and BO wrote the paper. All authors read and approved the final manuscript.

## Author details

<sup>1</sup> Mimetis BV, J.H. Oortweg 19, 2333 CH Leiden, The Netherlands. <sup>2</sup> Department of Cell and Chemical Biology, Leiden University Medical Centre, Einthovenweg 20, 2333 ZC Leiden, The Netherlands. <sup>3</sup> Biogen, 225 Binney Street, Cambridge, MA 02142, USA. <sup>4</sup> Yamaguchi University Graduate School of Medicine, Minamikogushi, Ube, Yamaguchi 7558505, Japan.

## Acknowledgements

The authors thank Dr. Thomas Cameron for his scientific input and Fang Qian for his help with reagent preparation and advice on MSD analysis.

## Competing interests

This study was funded by Biogen. Authors T.G., B.S., G.M., and B.O. are employees and shareholders of Biogen. Authors N.W., D.K., K.W., R.V., H.L., S.T. and P.V. are employees of MIMETAS and S.T. and P.V. are shareholders of that same company. The OrganoPlate® is a registered trademark of MIMETAS BV.

## Availability of data and materials

The materials and datasets are available from the authors on reasonable request.

## Consent for publication

Not applicable.

## Ethics approval and consent to participate

Not applicable.

## Funding

This work was partly supported by the CoSTREAM consortium, which has received funding from the European Union's Horizon 2020 research and innovation program under grant agreement No 667375, and the ADAPTED consortium, which has received funding from the Innovative Medicines Initiative 2 Joint Undertaking under grant agreement No 115975. This Joint Undertaking receives support from the European Union's Horizon 2020 research and innovation program and the European Federation of Pharmaceutical Industries and Associations.

## Publisher's Note

Springer Nature remains neutral with regard to jurisdictional claims in published maps and institutional affiliations.

Received: 17 May 2018 Accepted: 30 July 2018

Published online: 31 August 2018

## References

- Abbott NJ, Patabendige AAK, Dolman DEM, Yusof SR, Begley DJ. Structure and function of the blood–brain barrier. *Neurobiol Dis*. 2010;37:13–25.
- Wolburg H, Lippoldt A. Tight junctions of the blood–brain barrier: development, composition and regulation. *Vasc Pharmacol*. 2002;38:323–37.
- Daneman R. The blood–brain barrier in health and disease. *Ann Neurol*. 2012;72:648–72.
- Obermeier B, Daneman R, Ransohoff RM. Development, maintenance and disruption of the blood–brain barrier. *Nat Med*. 2013;19:1584–96.
- Wevers NR, de Vries HE. Morphogens and blood–brain barrier function in health and disease. *Tissue Barriers*. 2016;4:e1090524.
- Pardridge WM. Drug transport across the blood–brain barrier. *J Cereb Blood Flow Metab*. 2012;32:1959–72.
- Pardridge WM. Drug and gene delivery to the brain: the vascular route. *Neuron*. 2002;36:555–8.
- Freskgård PO, Urich E. Antibody therapies in CNS diseases. *Neuropharmacology*. 2017;120:38–55.
- Niewoehner J, Bohrmann B, Collin L, et al. Increased brain penetration and potency of a therapeutic antibody using a monovalent molecular shuttle. *Neuron*. 2014;81:49–60.
- Webster CI, Caram-Salas N, Haqqani AS, et al. Brain penetration, target engagement, and disposition of the blood–brain barrier-crossing bispecific antibody antagonist of metabotropic glutamate receptor type 1. *FASEB J*. 2016;30:1927–40.
- Zuchero JY, Chen X, Bien-Ly N, et al. Discovery of novel blood–brain barrier targets to enhance brain uptake of therapeutic antibodies. *Neuron*. 2016;89:70–82.
- Yu YJ, Zhang Y, Kenrick M, et al. Boosting brain uptake of a therapeutic antibody by reducing its affinity for a transcytosis target. *Sci Transl Med*. 2011;3:84ra44.
- Couch JA, Yu YJ, Zhang Y, et al. Addressing safety liabilities of TfR bispecific antibodies that cross the blood–brain barrier. *Sci Transl Med*. 2013;5:183ra57.
- Haqqani AS, Delaney CE, Brunette E, et al. Endosomal trafficking regulates receptor-mediated transcytosis of antibodies across the blood brain barrier. *J Cereb Blood Flow Metab*. 2018;38:727–40.
- Thom G, Hatcher J, Hearn A, et al. Isolation of blood–brain barrier-crossing antibodies from a phage display library by competitive elution and their ability to penetrate the central nervous system. *MAbs*. 2018;10:304–14.
- Weber F, Bohrmann B, Niewoehner J, et al. Brain shuttle antibody for Alzheimer's disease with attenuated peripheral effector function due to an inverted binding mode. *Cell Rep*. 2018;22:149–62.
- Ben-Zvi A, Lacoste B, Kur E, et al. Mfsd2a is critical for the formation and function of the blood–brain barrier. *Nature*. 2014;509:507–11.
- Villaseñor R, Ozmen L, Messaddeq N, et al. Trafficking of endogenous immunoglobulins by endothelial cells at the blood–brain barrier. *Sci Rep*. 2016;6:25658.
- Abbott NJ, Hughes CC, Revest PA, Greenwood J. Development and characterisation of a rat brain capillary endothelial culture: towards an in vitro blood–brain barrier. *J Cell Sci*. 1992;103(Pt 1):23–37.
- Biegel D, Pachter JS. Growth of brain microvessel endothelial cells on collagen gels: applications to the study of blood–brain barrier physiology and CNS inflammation. *In Vitro Cell Dev Biol Anim*. 1994;30:581–8.
- Gardner TW, Lieth E, Khin SA, et al. Astrocytes increase barrier properties and ZO-1 expression in retinal vascular endothelial cells. *Invest Ophthalmol Vis Sci*. 1997;38:2423–7.
- He Y, Yao Y, Tsirka SE, Cao Y. Cell-culture models of the blood–brain barrier. *Stroke*. 2014;45:2514–26.
- Wolff A, Antfolk M, Brodin B, Tenje M. In vitro blood–brain barrier models—an overview of established models and new microfluidic approaches. *J Pharm Sci*. 2015;104:2727–46.
- van der Helm MW, van der Meer AD, Eijkel JCT, van den Berg A, Segerink LI. Microfluidic organ-on-chip technology for blood–brain barrier research. *Tissue Barriers*. 2016;4:e1142493.
- Reichel A, Begley DJ, Abbott NJ. An overview of in vitro techniques for blood–brain barrier studies. *Methods Mol Med*. 2003;89:307–24.
- van der Meer AD, Orlova VV, ten Dijke P, van den Berg A, Mummery CL. Three-dimensional co-cultures of human endothelial cells and embryonic stem cell-derived pericytes inside a microfluidic device. *Lab Chip*. 2013;13:3562.
- Prabhakarpandian B, Shen MC, Nichols JB, et al. SyM-BBB: a microfluidic blood brain barrier model. *Lab Chip*. 2013;13:1093.
- Herland A, Van Der Meer AD, FitzGerald EA, Park TE, Sleeboom JF, Ingber DE. Distinct contributions of astrocytes and pericytes to neuroinflammation identified in a 3D human blood–brain barrier on a chip. *PLoS ONE*. 2016;11:1–21.
- Bang S, Lee SR, Ko J, et al. A low permeability microfluidic blood–brain barrier platform with direct contact between perfusable vascular network and astrocytes. *Sci Rep*. 2017;7:1–10.
- Sano Y, Shimizu F, Abe M, et al. Establishment of a new conditionally immortalized human brain microvascular endothelial cell line retaining an in vivo blood–brain barrier function. *J Cell Physiol*. 2010;225:519–28.
- Sano Y, Kashiwamura Y, Abe M, et al. Stable human brain microvascular endothelial cell line retaining its barrier-specific nature independent of the passage number. *Clin Exp Neuroimmunol*. 2013;4:92–103.
- Shimizu F, Sano Y, Tominaga O, Maeda T, Abe M, Kanda T. Advanced glycation end-products disrupt the blood–brain barrier by stimulating the release of transforming growth factor- $\beta$  by pericytes and vascular endothelial growth factor and matrix metalloproteinase-2 by endothelial cells in vitro. *Neurobiol Aging*. 2013;34:1902–12.
- Haruki H, Sano Y, Shimizu F, et al. NMO sera down-regulate AQP4 in human astrocyte and induce cytotoxicity independent of complement. *J Neurol Sci*. 2013;331:136–44.
- Shimizu F, Sano Y, Abe M, et al. Peripheral nerve pericytes modify the blood-nerve barrier function and tight junctional molecules through the secretion of various soluble factors. *J Cell Physiol*. 2011;226:255–66.
- Takeshita Y, Obermeier B, Cotleur A, Sano Y, Kanda T, Ransohoff RM. An in vitro blood–brain barrier model combining shear stress and endothelial cell/astrocyte co-culture. *J Neurosci Methods*. 2014;232:165–72.



36. Spampinato SF, Obermeier B, Cotleur A, et al. Sphingosine 1 phosphate at the blood brain barrier: can the modulation of S1P receptor 1 influence the response of endothelial cells and astrocytes to inflammatory stimuli? *PLoS ONE*. 2015;10:e0133392.
37. Takeshita Y, Obermeier B, Cotleur AC, et al. Effects of neuromyelitis optica-IgG at the blood-brain barrier in vitro. *Neurol Neuroimmunol Neuroinflamm*. 2017;4:e311.
38. Shimizu F, Schaller KL, Owens GP, et al. Glucose-regulated protein 78 autoantibody associates with blood-brain barrier disruption in neuromyelitis optica. *Sci Transl Med*. 2017;9:eaai9111.
39. Schindelin J, Arganda-Carreras I, Frise E, et al. Fiji: an open-source platform for biological-image analysis. *Nat Methods*. 2012;9:676–82.
40. Trietsch SJ, Israëls GD, Joore J, Hankemeier T, Vulto P. Microfluidic titer plate for stratified 3D cell culture. *Lab Chip*. 2013;13:3548–54.
41. Wevers NR, van Vught R, Wilschut KJ, et al. High-throughput compound evaluation on 3D networks of neurons and glia in a microfluidic platform. *Sci Rep*. 2016;6:38856.
42. van Duinen V, van den Heuvel A, Trietsch SJ, et al. 96 perfusable blood vessels to study vascular permeability in vitro. *Sci Rep*. 2017;7:18071.
43. Tietz S, Engelhardt B. Brain barriers: crosstalk between complex tight junctions and adherens junctions. *J Cell Biol*. 2015;209:493–506.
44. Armstrong JK, Wenby RB, Meiselman HJ, Fisher TC. The hydrodynamic radii of macromolecules and their effect on red blood cell aggregation. *Biophys J*. 2004;87:4259–70.
45. Hawe A, Hulse WL, Jiskoot W, Forbes RT. Taylor dispersion analysis compared to dynamic light scattering for the size analysis of therapeutic peptides and proteins and their aggregates. *Pharm Res*. 2011;28:2302–10.
46. Trietsch SJ, Naumovska E, Kurek D, et al. Membrane-free culture and real-time barrier integrity assessment of perfused intestinal epithelium tubes. *Nat Commun*. 2017;8:262.
47. Sade H, Baumgartner C, Hugenmatter A, Moessner E, Freskgård P-O, Niewoehner J. A human blood-brain barrier transcytosis assay reveals antibody transcytosis influenced by pH-dependent receptor binding. *PLoS ONE*. 2014;9:e96340.
48. Goldbaum FA, Cauerhff A, Velikovskiy CA, Llera AS, Riottot MM, Poljak RJ. Lack of significant differences in association rates and affinities of antibodies from short-term and long-term responses to hen egg lysozyme. *J Immunol*. 1999;162:6040–5.
49. Gumbleton M, Audus KL. Progress and limitations in the use of in vitro cell cultures to serve as a permeability screen for the blood-brain barrier. *J Pharm Sci*. 2001;90:1681–98.
50. Cho C-F, Wolfe JM, Fadzen CM, et al. Blood-brain-barrier spheroids as an in vitro screening platform for brain-penetrating agents. *Nat Commun*. 2017;8:15623.
51. Lippmann ES, Azarin SM, Kay JE, et al. Derivation of blood-brain barrier endothelial cells from human pluripotent stem cells. *Nat Biotechnol*. 2012;30:783–91.
52. Lippmann ES, Al-Ahmad A, Azarin SM, Palecek SP, Shusta EV. A retinoic acid-enhanced, multicellular human blood-brain barrier model derived from stem cell sources. *Sci Rep*. 2014;4:4160.
53. Katt ME, Linville RM, Mayo LN, Xu ZS, Searson PC. Functional brain-specific microvessels from iPSC-derived human brain microvascular endothelial cells: the role of matrix composition on monolayer formation. *Fluids Barriers CNS*. 2018;15:7.
54. Hollmann EK, Bailey AK, Potharazu AV, Neely MD, Bowman AB, Lippmann ES. Accelerated differentiation of human induced pluripotent stem cells to blood-brain barrier endothelial cells. *Fluids Barriers CNS*. 2017;14:9.
55. Ribocco-Lutkiewicz M, Sodja C, Haukenfrers J, et al. A novel human induced pluripotent stem cell blood-brain barrier model: applicability to study antibody-triggered receptor-mediated transcytosis. *Sci Rep*. 2018;8:1873.
56. Biegel D, Pachter JS. Growth of brain microvessel endothelial cells on collagen gels: applications to the study of blood-brain barrier physiology and CNS inflammation. *In Vitro Cell Dev Biol Anim*. 1994;30A:581–8.
57. Hopkins AM, DeSimone E, Chwalek K, Kaplan DL. 3D in vitro modeling of the central nervous system. *Prog Neurobiol*. 2015;125:1–25.
58. Naik P, Cucullo L. In vitro blood-brain barrier models: current and perspective technologies. *J Pharm Sci*. 2012;101:1337–54.
59. Dewey CF, Bussolari SR, Gimbrone MA, Davies PF. The dynamic response of vascular endothelial cells to fluid shear stress. *J Biomech Eng*. 1981;103:177.
60. Cucullo L, Hossain M, Tierney W, Janigro D. A new dynamic in vitro modular capillaries-venules modular system: cerebrovascular physiology in a box. *BMC Neurosci*. 2013;14:18.
61. Seebach J, Dieterich P, Luo F, et al. Endothelial barrier function under laminar fluid shear stress. *Lab Invest*. 2000;80:1819–31.
62. Siddharthan V, Kim YV, Liu S, Kim KS. Human astrocytes/astrocyte-conditioned medium and shear stress enhance the barrier properties of human brain microvascular endothelial cells. *Brain Res*. 2007;1147:39–50.
63. Cucullo L, Hossain M, Puvanna V, Marchi N, Janigro D. The role of shear stress in blood-brain barrier endothelial physiology. *BMC Neurosci*. 2011;12:40.
64. Reinitz A, DeStefano J, Ye M, Wong AD, Searson PC. Human brain microvascular endothelial cells resist elongation due to shear stress. *Microvasc Res*. 2015;99:8–18.
65. DeStefano JG, Xu ZS, Williams AJ, Yimam N, Searson PC. Effect of shear stress on iPSC-derived human brain microvascular endothelial cells (dhBMECs). *Fluids Barriers CNS*. 2017;14:20.
66. Zhang H, Chiao M. Anti-fouling coatings of poly(dimethylsiloxane) devices for biological and biomedical applications. *J Med Biol Eng*. 2015;35:143–55.
67. Wong I, Ho CM. Surface molecular property modifications for poly(dimethylsiloxane) (PDMS) based microfluidic devices. *Microfluid Nanofluidics*. 2009;7:291–306.
68. Jefferies WA, Brandon MR, Hunt SV, Williams AF, Gatter KC, Mason DY. Transferrin receptor on endothelium of brain capillaries. *Nature*. 1984;312:162–3.
69. Johnsen KB, Moos T. Revisiting nanoparticle technology for blood-brain barrier transport: unfolding at the endothelial gate improves the fate of transferrin receptor-targeted liposomes. *J Control Release*. 2016;222:32–46.
70. Farrington GK, Caram-Salas N, Haqqani AS, et al. A novel platform for engineering blood-brain barrier-crossing bispecific biologics. *FASEB J*. 2014;28:4764–78.
71. Ober RJ, Radu CG, Ghetie V, Ward ES. Differences in promiscuity for antibody-FcRn interactions across species: implications for therapeutic antibodies. *Int Immunol*. 2001;13:1551–9.
72. Garg A, Balthasar JP. Investigation of the Influence of FcRn on the distribution of IgG to the brain. *AAPS J*. 2009;11:553–7.
73. Moos T, Morgan EH. Restricted transport of anti-transferrin receptor antibody (OX26) through the blood-brain barrier in the rat. *J Neurochem*. 2001;79:119–29.